
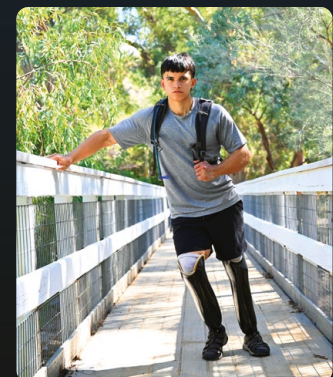
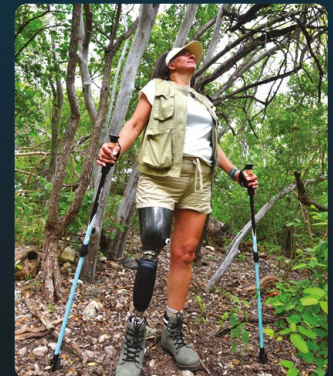


US Product Catalog

Components, Digital and Materials

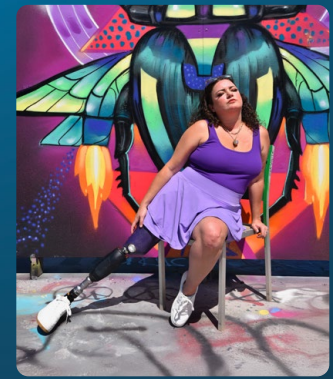

For over a century, PROTEOR has been driven by something deeper than technology—the belief that every person deserves to live fully, boldly, and without limits.

Our commitment to **#humanfirst** isn't just a motto; it's the heart of everything we do. We never stop pushing for what's possible, because we know every breakthrough means someone, somewhere, gets to walk a little farther, climb a little higher, or simply laugh a little louder.

We're proud to offer the most complete prosthetic portfolio in the industry, built not just on science, but on empathy and understanding. Our solutions are as unique as the people who wear them, designed so every prosthesis user can truly LIVE THE LIFE THEY LOVE.

We Are PROTEOR.

TABLE OF CONTENTS

L-Code Chart	4
SYNSYS	5
SYNSYS	7
Knees	10
PROTEOR QUATTRO 2™	13
Plié® 3	14
Symphony	15
HyTrek	16
Matik	17
1M Series	18
1M102 Series	19
1M112 Series	20
EasyRide	21
Ankles	22
PROTEOR Kinnex™ 2.0	25
Kinterra®	26
Feet	27
ShockWave®	31
RUSH ROGUE® 2	32
RUSH RAMPAGE®	33
PROTEOR Agilix™	34
DynAdapt	35
RUSH HiPro	36
Highlander®/ Highlander MAX	37
Sierra®	38
RUSH RAMPAGE® LP	39
ETHOS™ LP	40
Pacifica® LP	41
RUSH Chopart®	42
DynaStep	43
Gery	44
EasyRun	45
Hopper	46
RUSH H2O Collection	47
RUSH EVAQ8® Collection	48
EVAQ8 Collection	49
Foot Shells & Spectra Socks	50
Digital	52
Connectors & Components	59
Socket Technology	71
Materials	89
Ordering, Returns, Warranty Policies	104

**Easy online
ordering.
24 hours a day.
7 days a week.**

**Create your account at
shop.proteorusa.com**



L-CODE CHART

Feet	Item #	K-Level	L-Code(s)
Gery	1A200	K1-K2	L5972
PROTEOR Agilix	F15	K3	L5987, L5986
DynaStep® (w/multi-axis ankle®)	1A101	K3	L5981®, (L5986®)
DynAdapt	F10	K3	L5980, L5986
ETHOS LP	LP3	K3-K4	L5981, L5986
Highlander® & Highlander MAX	FS3	K3	L5981®
Pacifica LP	FS4	K3	L5981
RUSH Chopart®	CHO	K3-K4	L5976®
RUSH HiPro®	HiP	K3-K4	L5981®
RUSH RAMPAGE LP	RAMLP	K3-K4	L5981, L5986
RUSH RAMPAGE	RAM	K3-K4	L5987, L5986
RUSH ROGUE 2®	ROG2	K3-K4	L5987®, L5984®
ShockWave®	F20	K3-K4	L5987®, L5984®
Sierra	FS1	K3-K4	L5981, L5986
Vacuum Options			
RUSH HiPro EVAQ8	EVQH	K3-K4	L5981, L5781
RUSH RAMPAGE EVAQ8	EVRAM	K3-K4	L5987, L5986, L5781
RUSH RAMPAGE LP EVAQ8	EVQRAML	K3-K4	L5981, L5986, L5781
RUSH ROGUE 2 EVAQ8®	EVQR2	K3-K4	L5987®, L5984®, L5781
ShockWave EVAQ8®	F20-V2	K3-K4	L5987®, L5984®, L5781
PROTEOR Agilix EVAQ8	F15-V2	K3	L5987, L5986, L5781
Sierra EVAQ8	FS1-V2	K3-K4	L5981, L5986, L5781
Highlander EVAQ8 & Highlander MAX EVAQ8	FS3-V2	K3	L5981, L5781
DynAdapt EVAQ8	F10-V2	K3	L5980, L5986, L5781
ETHOS LP EVAQ8	LP3-V3	K3	L5981, L5986, L5781
Water Feet			
RUSH HiPro H2O	H2H	K3-K4	L5981
RUSH RAMPAGE H2O	H2RAM	K3-K4	L5987, L5986
RUSH RAMPAGE LP H2O	H2RAMLP	K3-K4	L5981, L5986
RUSH ROGUE 2 H2O®	H2R2	K3-K4	L5987®, L5984®
Sport Specific			
EasyRun	1A920	K4	L5980
Hopper	HLO	K4	L5999

Microprocessor Full-Leg System	Item #	K-Level	L-Code(s)
SYNSYS	1P700LEG	K2-K3	L5828®, L5856®, L5845®, L5848®, L5850®, L5925®, L5973®
Hydraulic Ankle			
Kinterra®	RM3	K3	L5981®, L5968®
Kinterra EVAQ8	RM3-V2	K3	L5981, L5968, L5781
Microprocessor Ankle			
PROTEOR Kinnex 2.0®	F14-N2	K3	L5973®
Microprocessor Knees			
Plié 3®	KS3/KD3	K2-K4	L5828®, L5856®, L5845®, L5848®, L5850®
PROTEOR QUATTRO 2®	QNX0010/ QNX0011	K3-K4	L5856®, L5828®, L5848®, L5845®, L5850®, L5925®
Mechanical Knees			
Single-Axis CF Locking	1M01/1M10	K1	L5811
4-Axis Adjustable Linkage	1M05	K1-K2	L5611, L5850
4-Axis Knee Long Linkage	1M102/HD	K1-K2	L5611, L5850
4-Axis Knee Long Linkage - w/Lock	1M102V/HD	K1-K2	L5611, L5850, L5925
Compact 4-Axis Knee	1M112/1M113	K4 for peds	L5611, L5850
Matik 4-Axis Pneumatic Knee	1P200/KD	K3	L5840, L5850
HyDeal II Hydraulic 4-Bar Knee - Grey Pyramid Top	1P120/KD	K3-K4	L5613, L5850
HyDeal II Hydraulic 4-Bar Knee - Orange Pyramid Top	1P121/KD	K3-K4	L5613, L5850
HyTrek	1P131/KD	K3-K4	L5828, L5845, L5850, L5925
Symphony - 6-Axis Hydraulic Knee w/Lock	NK-6/SH+L	K3	L5814, L5845, L5848, L5850, L5925
Symphony - 6-Axis Hydraulic Knee	NK-6/SH	K3	L5814, L5845, L5848, L5850
EasyRide Sport Knee	1P650	K4	L5830, L5845, L5848, L5850, L5925, L5930



The responsibility for accurate coding lies with the patient care facility that selects the product, fits the patient, and bills for the service. The product/device "Supplier", defined as an O&P practitioner, O&P patient care facility, or DME Supplier, assumes full responsibility for accurate billing of PROTEOR products. It is the Supplier's responsibility to determine medical necessity; ensure coverage criteria is met; and submit appropriate HCPCS codes, modifiers, and charges for services/products delivered. It is also recommended that Supplier's contact insurance payer(s) for coding and coverage guidance prior to submitting claims. PROTEOR Coding Suggestions and Reimbursement Guides are based on reasonable judgment and are not recommended to replace the Supplier's judgment. These recommendations may be subject to revision based on additional information or alpha-numeric system changes.



S4N/S4S

Standing is something most people never think about. The simple act of rising, balancing, and taking that first step forward feels automatic. Movement at this scale requires coordination, strength, and trust across the entire body.

Full-leg systems support alignment from hip to ground and distribute force across the entire structure. Each component works in coordination with the next, creating continuity through every phase of movement. The experience feels connected, stable, and whole.

SYNSYS

TABLE OF CONTENTS

Get more info by clicking on pictures



synsys
MICROPROCESSOR FULL-LEG SYSTEM




CONTACT US

ORDER



SYN SYS

MICROPROCESSOR FULL-LEG SYSTEM

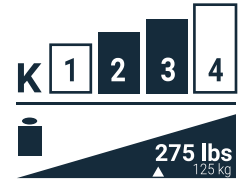

FULL-LEG SYSTEM

AN ENTIRELY NEW ANGLE

Only the **SYN SYS** Microprocessor Full-Leg System offers **triple flexion movements**, unlocking an entirely new angle on performance for K2 and low to moderate impact K3 ambulators.

An Entirely New Angle On:

Motion: Only **SYN SYS** provides hip-knee-ankle **triple flexion** with 125° knee flexion and 42° ankle ROM to move freely like never imagined.

Stability: **92% greater flat-foot time** for greater stability and **2x more toe clearance** for optimal security†. Intuitively walk downstairs with the full foot flat on the step ; and rely on its unique stumble recovery assistance.†

Toe Clearance: Active dorsiflexion during swing

Life: Up to 9 shoe options plus barefoot walking, 10+ days of battery life, and 13 protective cover colors to choose from.

Ordering Information:

- Certification is required before placing orders (for new patient trials and purchases).
- Certification courses are provided by the PROTEOR Clinical team and are free of charge.
- Patients are required to trial prior to purchasing.



L5828 L5856
L5845 L5848
L5973 L5850
L5925



reddot award
product design 2023



Programmable via
SYN SYS Soft App
for Prosthetists



Adjustable via
SYN SYS App
for Users



†Requena, C et al. (2023, September 18). Increased of stability and security in transfemoral amputees with a knee-ankle synergistic system. Presented at ESMAC. Athens, Greece.

DOWNLOAD PRODUCT SHEET



CONTACT US

ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com





SYNSYS

MICROPROCESSOR FULL-LEG SYSTEM

K

1

2

3

4

i

275 lbs

125 kg

FULL-LEG SYSTEM



- L5828L5845L5973L5925
- L5856L5848L5850



reddot award
product design 2023



TECHNICAL FEATURES

Item N°	1P700LEG
Proximal Connector	Pyramid
Weight	3220g / 7.1lbs
Knee-to-ground Height	16.9" to 21.7" / 430 to 550mm
Knee Flexion	125°
Ankle ROM	42° (20° dorsiflexion / 22° plantarflexion)
Max Patient Weight	275lbs / 125kg
Activity Level	K2 and low to moderate impact K3
Battery Life	>10 days depending on activity. Complete recharge time of 6 hours
CP Programming	iPad Pro 12.9" (5th generation), Windows 10 and 11
Patient App	iOS and Android
Water Resistance	IP54
Foot Sizes	23cm to 30cm
Foot Shell Colors	Light, Medium, Black

*Not suitable for bilateral transfemoral

Table of module availability and selection by patient's weight and foot size:

PATIENT WEIGHT	FOOT MODULE	FOOT SIZES							
		23	24	25	26	27	28	29	30
45-59 kg 99-130 lbs	1								
60-74 kg 131-163 lbs	2								
75-89 kg 164-196 lbs	3								
90-104 kg 197-229 lbs	4								
105-125 kg 230-275 lbs	5								

SYNSYS Protection Covers

1P700218-01-S
METALLIC GREY

1P700218-02-S
CREAMY WHITE

1P700218-03-S
MATTE BLACK

1P700218-06-S
MEDIUM GREY

1P700218-08-S
LIGHT FLESH

1P700218-09-S
MEDIUM FLESH

1P700218-10-S
DARK FLESH

1P700218-13-S
RUBY RED

1P700218-16-S
CAPRI BLUE

1P700218-17-S
METALLIC BLUE

1P700218-19-S
CAMELEON

1P700218-20-S
METALLIC VIOLET

1P700218-23-S
RAW MAT. COLOR

Lightweight, hygienic, and easy-to-maintain protection that is included with each purchase.

Available in **12 colors**, plus a raw material finish for patient customization.



DOWNLOAD PRODUCT SHEET



CONTACT US

ORDER



Knees

Every step forward is a decision. A willingness to bend without breaking, to rise again, to keep moving even when the path asks more of us.

Knees are designed to absorb force, provide stability, and support balance when direction changes. They make it possible to climb, kneel, turn, and stand tall again. With the right support, movement feels controlled, strong, and steady.

KNEES

TABLE OF CONTENTS



PROTEOR
QUATTRO 2
MICROPROCESSOR KNEE



PLIÉ3
MICROPROCESSOR KNEE



Symphony



HYTREK



matik



1M
SERIES



1M102
SERIES



1M112
SERIES



EASY
RIDE

[CONTACT US](#)

[ORDER](#)

[f](#) [i](#) [v](#) [in](#) **#HUMANFIRST**

© 2026 by PROTEOR | All rights reserved.

us.proteor.com






PROTEOR

QUATTRO 2™

MICROPROCESSOR KNEE

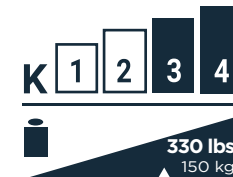

KNEES



L5828 L5850
L5845 L5856
L5848 L5925

EXPERIENCE THE EVOLUTION.

QUATTRO 2 doesn't just sense your movement. It adapts in real time. Every step feels like your own.

COMFORT IN MOTION

Traditional MPKs rely on dual-valve systems that open and close between steps, creating lag. QUATTRO 2 runs on a single bidirectional valve: continuous resistance, twice the speed, zero lag. Motion feels natural. Comfort feels automatic.

FASTER THAN YOUR BIOLOGICAL REFLEXES

A processor 10x faster than the original QUATTRO. Three high-precision sensors reading your gait in real time. Every transition, instinctive.

- Angle sensor 10x more sensitive than traditional MPKs
- Sampling rate 2x greater than conventional designs

ALWAYS ONE STEP AHEAD

Your knee keeps getting smarter. Remote firmware updates, powered by real-world usage and Snapshot software, unlock new features and sharpen performance over time. (Patients opt in to share data.)

BUILT AROUND YOU

Engineered to fit your life, not the other way around.

- The most compact MPK on the market: just 8.5"
- Now rated for 330 lbs / 150 kg: a stronger frame, built for more patients
- Adaptive gait tuning, right from your phone

TECHNICAL FEATURES

Item N°.	QNX0010	QNX0011
Proximal Connector	Pyramid	Threaded Top
Weight	1655g / 3.65lbs	
Build Height	21.6 cm / 8.5 in	
Knee Flexion	135°	
Maximum Patient Weight	330lbs / 150kg	
Activity Level	K3 / K4	
Battery Life	2-3 Days	
Water Resistance	IP68	
Booster Adapter	ACC0010	
Warranty	36 Months	

WHAT'S NEW WITH 2

Stair Ascent | IP68 Rating | 330 LB Weight Limit
2-3 Day Battery Life | 10x Faster Processor | USB-C

Programmable via
GaitLab App
for Prosthetists
(iOS and Android)



Adjustable via
Freedom Innovations App
for Users
(iOS and Android)



DOWNLOAD PRODUCT SHEET



CONTACT US

ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com



A Fundamental Fit

Empowering patients with dependable performance and essential functionality. The **Plié 3**® Microprocessor Knee is a dependable knee that’s versatile enough for anywhere your K2 and low, mod, high K3 and K4.

- Streamlined design with rugged internal components and modern aesthetics
- Quickly responds, adapting to a variety of gait speeds and changes in user activities
- Customizable swing release parameters
- Manual adjustments available to optimize activities
- Interchangeable batteries for on-the-go convenience
- IP67 rating for full submersion in water
- Generation of cadence reports



L5828L5850

L5845L5856


L5848

TECHNICAL FEATURES

Item N°.	KS3-00-KNEE1-KT	KD3-00-KNEE2-KT
Proximal Connector	Pyramid	Threaded Top
Weight	1235g / 2.7lbs	1243g / 2.4lbs
Build Height	9.25" / 235mm	8.75" / 223mm
Knee Flexion	125°	117°
Maximum Patient Weight	275lbs / 125kg moderate activity 220lbs / 100kg high activity	
Activity Level	K2 / K3 / K4	
Battery Life	24 hours per battery (2 interchangeable batteries included)	
Water Resistance	IP67	

Fairing Options



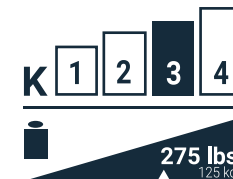

KNEES

Fluidity, Functionality, Safety; A Perfect Composition


The **Symphony** Knee is designed with a unique, 6-Bar polycentric structure that provides an extremely high level of stability during stance phase. The innovative p-MRS system, automatically detects gait status and controls the stability of the knee joint accordingly. The hydraulic cylinder provides a smooth swing phase and works with the **Symphony's** adjustable stance flexion feature to enhance the patient's ambulatory experience and provide a comfortable, confident gait to the user.

- 6-Bar hydraulic knee
- Adjustable stance flexion
- Only polycentric knee with an optional manual lock
- High level of stability due to geometric locking of the knee through toe off
- Low build height
- Unique hydraulic cylinder design allows easy flexion initiation



TECHNICAL FEATURES

Item N°.	NK-6	NK-6+L	NK-6 SH	NK-6 SH+L
Proximal Connection	Pyramid		Threaded Top	
Weight	920g / 2.0lbs	970g / 2.1lbs	960g / 2.1lbs	1010g / 2.2lbs
Build Height	7.7" / 197mm		7.5" / 191mm	
Knee Flexion	170°			
Maximum Patient Weight	275lbs / 125kg 220lbs / 100kg for Hip Prosthesis & Highly Active Use			
Activity Level	K3			
Selective Lock	No	Yes	No	Yes

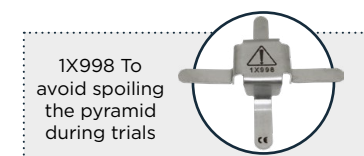[DOWNLOAD PRODUCT SHEET](#)[CONTACT US](#)[ORDER](#)

[f](#) [i](#) [v](#) [in](#) #HUMANFIRST


© 2026 by PROTEOR | All rights reserved.

us.proteor.com





HYTREK

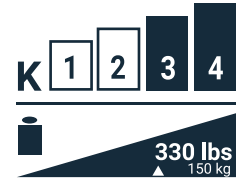

KNEES

Single Axis With Hydraulic Regulation of Stance and Swing Phases And Manual Lock

The **Hytrek** is a durable knee for dynamic transfemoral amputees. Navigate everyday life activities on all types of terrain; including slopes and stairs. Offers smooth ambulation and the ability to safely descend stairs step over step.


- Adjustable hydraulic resistance for supported stair decent
- Optional manual lock
- Extension assist
- Requires 34mm Pylon Distal Connectivity (sold separately)



- STANCE:**
Adjustment of sensitivity for switching to braking mode.
- RESISTANCE:**
Adjustment of braking power.
- LOCK:**
Manual hydraulic lock switch.
- FLEXION:**
Adjustment of flexion angle.
- EXTENSION:**
Dynamic control of swing phase and terminal impact.

TECHNICAL FEATURES

Item N°.	1P131	1P131-KD
Proximal Connector	Pyramid	Threaded Top
Weight	1265g / 2.8lbs	1285g / 2.8lbs
Build Height	7.0" / 177mm	6.9" / 174mm
Knee Flexion	120°	
Distal Connector	34mm pylon	
Maximum Patient Weight	330lbs / 150kg	
Activity level	K3 / K4	
Water Resistance	Weather resistant	
Materials	Aluminum, Stainless steel, Titanium	

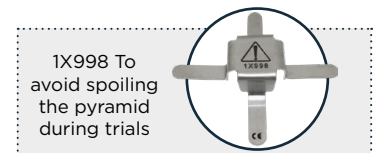[DOWNLOAD PRODUCT SHEET](#)[CONTACT US](#)[ORDER](#)

[f](#) [i](#) [v](#) [in](#) #HUMANFIRST


© 2026 by PROTEOR | All rights reserved.

us.proteor.com





matik

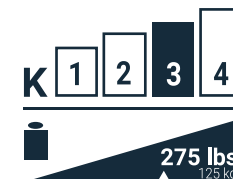

KNEES

An Affordable 4-axis Pneumatic Knee For Patients Up To 275lbs


Simply and effectively designed to provide high function with low maintenance. The easy swing of the pneumatic unit and the reduced weight of the knee give the patient a heightened sensation of fluidity and comfort. To support safety during swing phase and reduce the risk of a fall, the toe clearance during swing reaches up to 25mm, one of the largest clearances for this type of knee. **Matik** has both flexion and extension adjustments which address terminal impact and heel rise.

- Pneumatic knee with an adjustable extension assist
- Secure 4-axis geometry
- Safety and security through generous toe clearance during swing phase
- Compact and lightweight



TECHNICAL FEATURES

Item N°.	1P200	1P200-KD
Proximal Connector	Pyramid	Threaded Top
Weight	743g / 1.6lbs	
Build Height	6.9" / 175mm	6.9" / 175mm
Knee Flexion	170°	
Distal Connector	Pyramid	
Maximum Patient Weight	275lbs / 125kg	
Activity Level	Low to moderate K3	
Materials	Aluminum, Stainless steel	

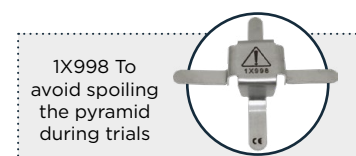


CONTACT US

ORDER



1M SERIES

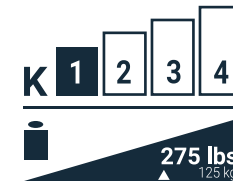

KNEES

Convenient Lock Ensures Patient Safety While Walking


- Joint friction adjustment
- Lock with adjustable stop
- Socket rotation adjustment: $\pm 15^\circ$
- 1M01 / 1M10 require 30mm pylon distal connectivity
- 1M01-HD / 1M10-HD require 34mm pylon distal connectivity




TECHNICAL FEATURES

Item N°.	Max. Patient Weight	Weight	TH	BH	H	Materials	Release lever	Flexion	Distal Connector
1M01	220lbs/ 100kg	345g/ 0.8lbs	3.9"/ 99mm	2.3"/ 58mm	0.5"/ 13mm	Aluminum Steel	No	130°	Ø 30mm
1M10		355g/ 0.8lbs					Yes		
1M01-HD	275lbs/ 125kg	375g/ 0.8lbs	4.0"/ 102mm	2.3"/ 59mm	0.6"/ 14mm		No		Ø 34mm
1M10-HD							Yes		

1M01 / 1M10

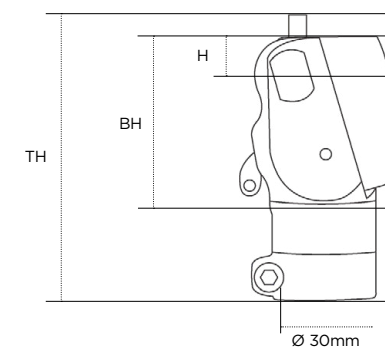

1M01-HD / 1M10-HD

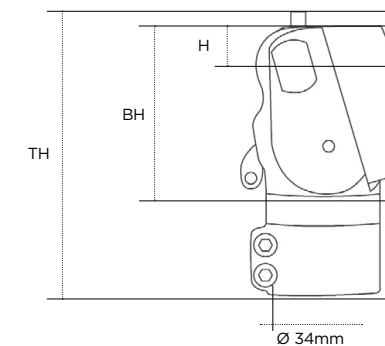

PROXIMAL ADAPTERS REQUIRED TO COMPLETE YOUR ORDER



1K40



1K160-HD



1K160



1K163



CONTACT US

ORDER

f i y t #HUMANFIRST


© 2026 by PROTEOR | All rights reserved.

us.proteor.com





1M102 SERIES

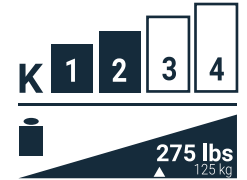

KNEES



Unique Alignment and Axis Geometry Offers Increased Safety During Stance Phase

- Socket rotation adjustment: $\pm 15^\circ$
- Friction adjustment
- Extension assist adjustment
- Optional lock available
- 1M02 / 1M102V require 30mm pylon distal connectivity
- 1M02-HD / 1M102V-HD require 34mm pylon distal connectivity



TECHNICAL FEATURES

Item N°.	Max. Patient Weight	Weight	TH	BH	H	Material	Lock	Flexion	Distal Connector
1M102	220lbs/ 100kg	565g/ 1.2lbs	7.2"/ 184mm	4.4"/ 113mm	0.5"/ 14mm	Aluminum	No	160°	Ø 30mm
1M102V		575g/ 1.3lbs					Yes		
1M102-HD	275lbs/ 125kg	500g/ 1.1lbs					No		Ø 34mm
1M102V-HD		525g/ 1.2lbs					Yes		
1X110	(required for 1M102V / 1M102V-HD) Socket Mounted Mechanical Knee Lock Release								

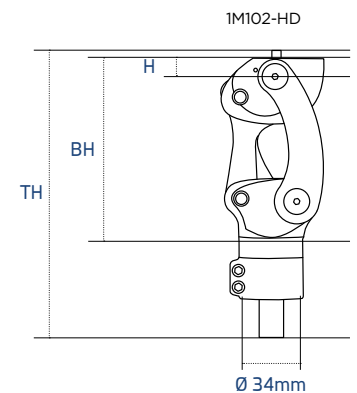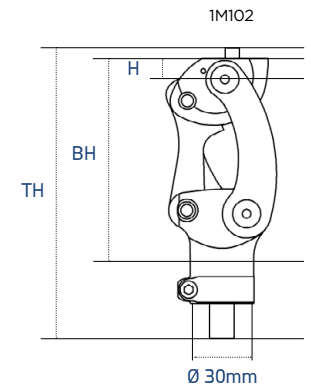

1X110

PROXIMAL ADAPTERS REQUIRED TO COMPLETE YOUR ORDER



1K40



1K160-HD



1K160



1K163



CONTACT US

ORDER



1M112 SERIES

K

1

2

3


4

176 lbs
80 kg

KNEES

Pediatric Compact 4-axis Knee

- Socket rotation adjustment: $\pm 15^\circ$
- Friction adjustment
- Extension assist adjustment
- Safety in stance phase ensured by alignment and axis geometry
- The 4-axis and linkage geometry provides a shortening of the leg segment during swing phase
- 1M112 requires 22mm pylon distal connectivity
- 1M113 requires 30mm pylon distal connectivity

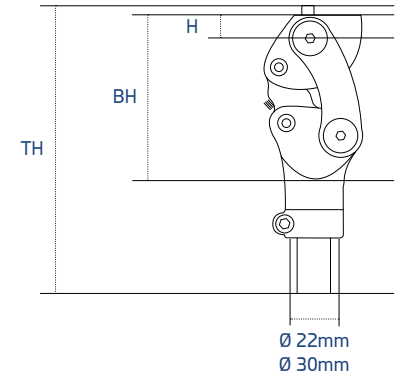

TECHNICAL FEATURES


	Item N°.	Max. Patient Weight	Weight	TH	BH	H	Material	Flexion	Distal Connector
Peds	1M112	99lbs / 45kg	320g / 0.7 lbs	6.1" / 155mm	3.0" / 75mm	0.4" / 11mm	Aluminum	165°	Ø 22mm
	1M113	176lbs / 80kg							Ø 30mm

PROXIMAL ADAPTERS REQUIRED TO COMPLETE YOUR ORDER



[CONTACT US](#)

[ORDER](#)

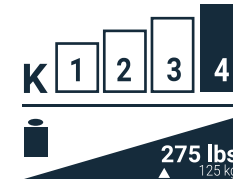

KNEES

High Impact Multi-sport Knee

The **EasyRide** knee was created in collaboration with adaptive athletes from around the world to ensure optimal performance capabilities for extreme, high impact, and all-terrain sports such as surfing, snowboarding, skateboarding, cycling, skiing, and more. The **EasyRide** lets you get back to the activities you LOVE without limitations.

Durability and Safety

- The FOX® shock absorber smoothly reproduces the natural motion of the knee and ensures smooth energy return
- Aeronautical-grade aluminum alloy and custom hardening treatment withstands most extreme conditions including saltwater and cold temperatures
- Lightweight, sturdy design

Adaptability

- Customizable for optimal function depending on type of sport, requiring knee flexion
- Adjusts to 12° pre-flexion / 60° max; biking position adjusted with 0° pre-flexion / 130° flexion max allowing for both good load transfer to the pedal and enough flexion for pedaling



TECHNICAL FEATURES

Item N°.	1P650
Proximal Connector	Pyramid
Distal Connector	Pyramid
Activity level	K4
Weight	1264g / 2.3lbs
Build Height / Full Height Position A	9" - 10" / 227mm - 257mm
Position A Pre-flexion / Flexion / Pression	12° / 60° max / 4 to 20 bars
Build Height / Full Height Position B	9" - 10" / 230mm - 260mm
Position B Pre-flexion / Flexion / Pression	0° / 130° max / 1 to 6 bars
Maximum Patient Weight	275lbs / 125kg

- Paddle-boarding
- Surfing
- Wind Surfing
- Water Skiing
- Wake Surfing
- Skiing (Downhill and Cross Country)
- Snowboarding
- Skateboarding
- Mountain Biking
- Road Biking
- And more...



DOWNLOAD PRODUCT SHEET



CONTACT US

ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com





Ankles

Every step we take is an act of trust. Trust that the ground will meet us, trust that our body will carry us forward, trust that movement **still belongs to us.**

Ankles absorb impact, adapt to uneven ground, and steady the body when life doesn't move in straight lines. With responsive support and natural articulation, each step feels controlled, fluid, and confident — allowing movement to continue without hesitation.

ANKLES

TABLE OF CONTENTS

Get more info by clicking on pictures



PROTEOR
KINNEX™ 2.0
MICROPROCESSOR ANKLE/FOOT SYSTEM



KINTERRA®



CONTACT US

ORDER


Designed for my real life... without compromise

The **PROTEOR Kinnex 2.0** integrates the world’s fastest responding microprocessor ankle/foot technology and carbon fiber footplate to provide low to moderate impact K3 ambulators the stability they expect — instantly. Unmatched ground compliance, improved comfort, and heel height adjustments keep users firmly grounded and connected, wet or dry.

- IP67 rated for water submersion up to 1m for 30 mins
- Sophisticated microprocessor controlled resistance helps optimize stability whenever changes in terrain or cadence occurs
- Accommodates heel heights from 0-2" and allows user to program up to 100 shoes
- 30° of microprocessor-controlled ankle range of motion: 20° plantarflexion / 10° dorsiflexion allows the ankle and foot to be ideally positioned in relation to the ground
- Split keel and heel for medial/lateral terrain adaptability
- Booster battery pack available for on-the-go charging



PDAC
VERIFIED
L5973

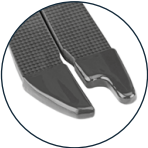

Sandal-toe available for sizes 24-28cm.

Programmable via
GaitLab App
for Prosthetists
(iOS and Android)



Adjustable via
Freedom Innovations App for Users
(iOS and Android)



TECHNICAL FEATURES

Item N°.		F14-N2
Weight (26 Cat 4, module with foot shell and spectra sock)		1641g / 3.6lbs
Build Height*	24 - 25cm	7.2" / 184mm
	26 - 28cm	7.3" / 186mm
	29 - 30cm	7.3" / 186mm
Standard Foot Shell (No CAP)		FTC-3K-1
Sandal Toe Foot Shell (No CAP)		FTC-2K-1

Spectra Sock and Heel Bumper Kit included.

Foot Shell must be ordered separately.

*Build height based on size 24cm, 26cm, 29cm Cat 4 foot module with foot shell, spectra sock, and 10mm heel height. Sandal-toe available for sizes 24-28cm.

Module selection according to activity and weight								
Weight/ Activity level	100-115lbs 44-52kg	116-130lbs 53-59kg	131-150lbs 60-68kg	151-170lbs 69-77kg	171-195lbs 78-88kg	196-220lbs 89-100kg	221-255lbs 101-116kg	256-275lbs 117-125kg
Low	1	1	2	3	4	5	6	7
Moderate	1	2	3	4	5	6	-	-

DOWNLOAD PRODUCT SHEET




CONTACT US

ORDER



KINTERRA®

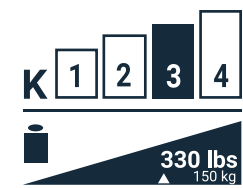

ANKLES

Now Available in
EVAQ8®



PDAC
VERIFIED
L5981 & L5968

Available in Sandal
Toe only

Redefining the Ankle Experience

The time-tested hydraulic ankle that continues to raise the bar for quality, stability, range of motion, and comfort.

- 17° Range of motion (2° dorsiflexion / 15° plantarflexion) for terrain adaptation and sitting comfort
- 330lb (150kg) user weight limit
- Submersible in fresh water up to 1m for 30 mins
- Completely redesigned hydraulic ankle and EnduraCore® full length footplate with split keel and heel
- Available in Sandal Toe only

TECHNICAL FEATURES

Item N°.		RM3
Weight (26 Cat 4, module with foot shell and spectra sock)		794g / 1.75lbs
Build Height*	22 - 25cm	4.65" / 118mm
	26 - 28cm	4.67" / 119mm
	29 - 30cm	4.81" / 122mm
Heel Height		3/8" / 10mm
Sandal Toe Foot Shell (No CAP)		FTC-2K-1

Black Spectra sock and
Heel Bumper Kit included.
Foot Shell must be
ordered separately.

*Build height based on size 23cm, 26cm, 29cm Cat 4 foot module with foot shell, spectra sock, and 10mm heel height.

Category selection according to activity and weight									
Weight/ Activity level	100-115lbs 44-52kg	116-130lbs 53-59kg	131-150lbs 60-68kg	151-170lbs 69-77kg	171-195lbs 78-88kg	196-220lbs 89-100kg	221-255lbs 101-116kg	256-290lbs 117-132kg	291-330lbs 133-150kg
Low	1	1	2	3	4	5	6	7	8
Moderate	1	2	3	4	5	6	7	8	-

The Kinterra 2.0 has been discontinued. Please see pg. 50 for replacement footshells for Kinterra 2.0

DOWNLOAD PRODUCT SHEET

€

CONTACT US

ORDER



Feet

Movement begins at the ground. With contact and balance. Every shift in weight is a conversation between body and surface, a language that keeps us steady and moving.

Feet adapt, stabilize, and respond in real time. Impact is cushioned, terrain is negotiated, and each shift in direction feels supported. The result is movement that feels natural, secure, and entirely our own.

FEET

PRODUCT LINKS

Get more info by clicking on pictures



SHOCKWAVE®



RUSH
ROGUE



RUSH
RAMPAGE



AGILIX™



DYNADAPT



RUSH
HiPro



HIGHLANDER®
HIGHLANDER MAX



SIERRA®



RUSH
RAMPAGE



ETHOS™ LP

NEW



PACIFICA™ LP

CONTACT US

ORDER

FEET

PRODUCT LINKS



FOOT SHELLS AND SPECTRA SOCKS

CONTACT US

ORDER

[f](#) [i](#) [v](#) [i](#) [n](#) #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com





SHOCKWAVE®

K

1

2

3

4

i


325 lbs
147 kg

FEET

Now Available in
EVAQ8®




PDAC
VERIFIED
L5987 & L5984

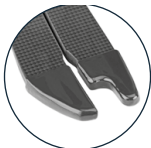

Sandal-toe available
for sizes 23-28cm.

Pure Excitement

Confidence for an active lifestyle—and the freedom to move as responsively as K3/K4 ambulators want to—are now possible with the premier **ShockWave®** carbon fiber foot from PROTEOR. Provides a shockingly smooth experience during high-impact activities, due to an advanced, multi-functional Spherical Shock Unit for ultimate vertical shock absorption.

- **Move More Naturally:** Allows a more natural gait by providing ultimate shock absorption at heel strike, which eases stress on the residual limb.
- **Seamless Movement:** Axial rotation of +/- 15° reduces socket pressure and allows for quick pivoting.
- **Trust Every Step:** Adds confidence when navigating diverse terrains with a split keel and heel providing 18° of inversion/eversion - improving stability, movement fluidity and comfort.
- **Lasting Reliability:** Lighter, more reliable foot due to a bolt-less footplate design that eliminates dead spots and weak structural areas.
- **Enjoy More Activities:** Proprietary EnduraCore technology provides ideal strength, flexibility, and durability for high energy return with less fatigue.
- **Go More Places:** Available in sandal toe.
- **Water Resistant**

TECHNICAL FEATURES

Item N°.		F20-S3
Weight (26 Cat4, module with foot shell and spectra sock)		978g / 2.2lbs
Build Height*	23 - 25cm	6.9" / 176mm
	26 - 28cm	7.4" / 189mm
	29 - 30cm	7.8" / 200mm
Heel Height		3/8" / 10mm
Standard Foot Shell (No CAP)		FTC-3F-1
Sandal Toe Foot Shell (No CAP)		FTC-2F-1

Black Spectra Sock and Heel Bumper Kit included in each Foot Module delivery.

Foot Shell must be ordered separately.

*Build height based on size 23cm, 26cm, 29cm Cat 4 foot module with foot shell, spectra sock, and 10mm heel height. Sandal-toe available for sizes 23-28cm.

Module selection according to activity and weight									
Weight/ Activity level	100-115lbs 44-52kg	116-130lbs 53-59kg	131-150lbs 60-68kg	151-170lbs 69-77kg	171-195lbs 78-88kg	196-220lbs 89-100kg	221-255lbs 101-116kg	256-285lbs 117-130kg	286-325lbs 131-147kg
Low	1	1	2	3	4	5	6	7	8
Moderate	1	2	3	4	5	6	7	8	-
High	2	3	4	5	6	7	8	-	-

1X998 To
avoid spoiling
the pyramid
during trials




DOWNLOAD PRODUCT SHEET



CONTACT US

ORDER

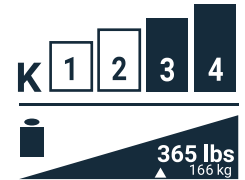

FEET

Now Available in
EVAQ8® and H2O



L5987 & L5984



NOW AVAILABLE!
Sandal Toe
Sizes 22-29cm

Pure Confidence

Steady and ready—with the ability to maintain durability during slow-load activities—are just some of the many benefits of **RUSH ROGUE® 2**. Remarkably reliable, this premier, fiberglass vertical shock foot from PROTEOR delivers increased dorsiflexion strength and the utmost flexibility.

- **Walk More Easily and Naturally:** Ideal dynamics design with a Heel Shock Bumper for maximum energy return at heel strike and 27° dorsi-plantar flexion amplitude.
- **Move More Ways:** Provides optimal durability to perform the way users want it to.
- **Seamless Movement:** Axial rotation of +/- 15° reduces socket pressure and allows for quick pivoting.
- **Go More Places:** Now available in sandal toe.
- **Waterproof** (fresh or saltwater)

TECHNICAL FEATURES


Item N°.		ROG2
Weight	(26 Cat 4, with foot shell and spectra sock)	1039g / 2.3lbs
Build Height*	22 - 24cm	6.1" / 155mm
	25 - 27cm	6.6" / 168mm
	28 - 29cm	6.7" / 171mm
Heel Height		3/8" / 10mm
Standard Foot Shell (With CAP)		FS
Sandal Toe Foot Shell (No CAP)		FTC-4R-1

*Build height based on size 23cm, 26cm, or 29cm Cat 4 foot module with spectra sock, foot shell and 10mm heel height.

Black Spectra Sock and Heel Wedge included in each Foot Module delivery.
Foot Shell must be ordered separately.

Module selection according to activity and weight									
Weight/ Activity level	0-105lbs 0-48kg	106-140lbs 49-64kg	141-175lbs 65-79kg	176-210lbs 80-95kg	211-245lbs 96-111kg	246-280lbs 112-127kg	281-315lbs 128-143kg	316-350lbs 144-159kg	351-365lbs 160-166kg
Low	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9
High	2	3	4	5	6	7	8	9	-

NOTE: Increasing category numbers indicate increased keel response. FOR BILATERAL PATIENTS: Increase by one category for added stability.

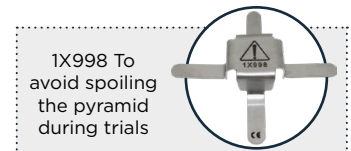

DOWNLOAD PRODUCT SHEET



CONTACT US


ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com



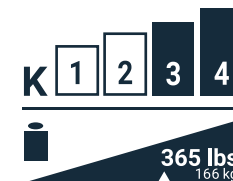

FEET

Now Available in
EVAQ8® and H2O



NOW AVAILABLE!
Sandal Toe
Sizes 22-29cm

Unbridled Capabilities

The smooth roll-through, increased comfort, and all-terrain capability of its predecessors, the **RUSH RAMPAGE®** boasts a more cosmesis-friendly option, with a lower overall product weight. This foot is the perfect solution for individuals seeking a blend of both cosmetic and high-performance capabilities experienced in the RUSH Foot Collection. Are you ready to RAMPAGE?

- Cosmesis-ready to meet your aesthetic preferences
- 95% energy return efficiency
- Smooth rollover with unparalleled energy return
- Extreme reliability for any use in any environment
- Waterproof (fresh and saltwater)

TECHNICAL FEATURES


Item N°.		RAM
Weight (26 Cat 4, with foot shell and spectra sock)		850g / 1.9lbs
Build Height*	22 - 24cm	6.6" / 168mm
	25 - 27cm	7.1" / 181mm
	28 - 29cm	7.2" / 184mm
Heel Height		3/8" / 10mm
Standard Foot Shell (With CAP)		FS
Sandal Toe Foot Shell (No CAP)		FTC-4R-1

*Build height based on size 23cm, 26cm, or 29cm Cat 4 foot module with spectra sock, foot shell and 10mm heel height.

Black Spectra Sock and Wedge included in each Foot Module delivery.
Foot Shell must be ordered separately.

Module selection according to activity and weight									
Weight/ Activity level	0-105lbs 0-48kg	106-140lbs 49-64kg	141-175lbs 65-79kg	176-210lbs 80-95kg	211-245lbs 96-111kg	246-280lbs 112-127kg	281-315lbs 128-143kg	316-350lbs 144-159kg	351-365lbs 160-166kg
Low	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9
High	2	3	4	5	6	7	8	9	-

NOTE: Increasing category numbers indicate increased keel response. FOR BILATERAL PATIENTS: Increase by one category for added stability.

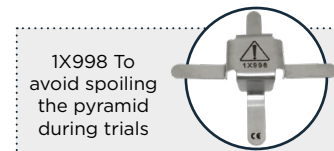

DOWNLOAD PRODUCT SHEET



CONTACT US

ORDER

f i y t #HUMANFIRST


© 2026 by PROTEOR | All rights reserved.

us.proteor.com






AGILIX™

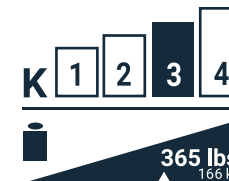

FEET

Now Available in
EVAQ8®

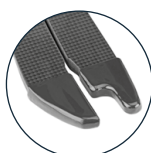

Sandal-toe available
for sizes 22-28cm.

Maximum Shock Value

The PROTEOR **Agilix™** is a multi-axial, shock absorbing, flexible foot system designed to manage loading impacts, reduce socket shear forces, and improve comfort while walking on any terrain. The ultra-lightweight design offers K3 ambulators the most exciting and comfortable daily experience.

- Water resistant
- Vertical shock absorption and carbon fiber flexibility optimized for low to moderate impact
- K3 ambulators
- Multi-axial function increases ground compliance and stability
- Streamlined, split keel and heel design results in improved long-term durability
- Full length dynamic heel attached at the toe ensures a seamless rollover, for superior comfort and a more symmetrical gait
- Lower profile and lighter weight than the competitive alternative, allowing users to do more and feel less tired at the end of the day

TECHNICAL FEATURES

Item N°.		F15
Weight (26 Cat 4, module with foot shell and spectra sock)		605g / 1.3lbs
Build Height*	22 - 25cm	5.5" / 140mm
	26 - 28cm	6.2" / 159mm
	29 - 31cm	6.8" / 173mm
Heel Height		3/8" / 10mm
Standard Foot Shell (No CAP)		FTC-3M-1
Standard Foot Shell (With CAP)		FTC-3M-0
Sandal Toe Foot Shell (No CAP)		FTC-2M-1

Black Spectra Sock and
Heel Bumper Kit included in
each Foot Module delivery.

Foot Shell must be
ordered separately.

*Build height based on size 23cm, 26cm, 29cm Cat 4 foot module with foot shell, spectra sock, and 10mm heel height. Sandal-toe available for sizes 22-28cm.

Module selection according to activity and weight										
Weight/ Activity level	100-115lbs 44-52kg	116-130lbs 53-59kg	131-150lbs 60-68kg	151-170lbs 69-77kg	171-195lbs 78-88kg	196-220lbs 89-100kg	221-255lbs 101-116kg	256-285lbs 117-130kg	286-325lbs 131-147kg	326-365lbs 148-166kg
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	-

1X998 To
avoid spoiling
the pyramid
during trials



CONTACT US


ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com



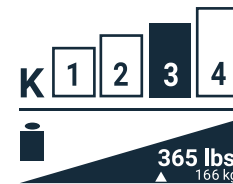

FEET

Now Available in
EVAQ8®

Fit, Finish, and Roll

The **DynAdapt** prosthetic foot is a slim profile carbon fiber foot system with a slender, anatomic design for easy fit and finish. Its multi-axial function provides maximum comfort and the uninterrupted strands of carbon fiber in the full-length keel and patients with effortless rollover and a more natural gait.


- Water-resistant
- A full length, unbolted sole plate ensures seamless stance phase roll-over for superior comfort
- Split keel and heel design increase ground compliance and improve stability
- Excellent dynamic response with 94% energy return efficiency

TECHNICAL FEATURES

Item N°.		F10
Weight (26 Cat 4, module with foot shell and spectra sock)		628g / 1.4lbs
Build Height*	22 - 25cm	6.1" / 155mm
	26 - 28cm	6.5" / 166mm
	29 - 31cm	7" / 178mm
Heel Height		3/8" / 10mm
Standard Foot Shell (No CAP)		FTC-3M-1
Standard Foot Shell (With CAP)		FTC-3M-0
Sandal Toe Foot Shell (No CAP)		FTC-2M-1

Black Spectra Sock and Heel Bumper Kit included in each Foot Module delivery.

Foot Shell must be ordered separately.

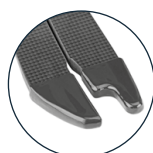

Sandal-toe available for sizes 22-28cm.

*Build height based on size 23cm, 26cm, 29cm Cat 4 foot module with foot shell, spectra sock, and 10mm heel height. Sandal-toe available for sizes 22-28cm.

Module selection according to activity and weight										
Weight/ Activity level	100-115lbs 44-52kg	116-130lbs 53-59kg	131-150lbs 60-68kg	151-170lbs 69-77kg	171-195lbs 78-88kg	196-220lbs 89-100kg	221-255lbs 101-116kg	256-285lbs 117-130kg	286-325lbs 131-147kg	326-365lbs 148-166kg
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	-
High	2	3	4	5	6	7	8	9	-	-

1X998 To avoid spoiling the pyramid during trials



CONTACT US

ORDER

Now Available in
EVAQ8[®] and H2O



NOW AVAILABLE!
Sandal Toe
Sizes 22-29cm

Virtually Indestructible

The **RUSH HiPro[®]** provides a realistic and responsive foot with smooth rollover. Made from an advanced engineered glass composite, wearers can feel confident in its durability with the added flexibility. The energy return of this foot enables amputees to easily traverse the most aggressive terrain while experiencing zero dead spots.

- Highly dynamic
- High energy return
- Smooth rollover
- All-terrain and waterproof (fresh or saltwater)

TECHNICAL FEATURES

Item N°.		HiP
Weight (26 Cat 4, with foot shell and spectra sock)		920g / 2.0lbs
Build Height*	22 - 24cm	6.5" / 165mm
	25 - 27cm	6.9" / 175mm
	28 - 29cm	7.1" / 181mm
Heel Height		3/8" / 10mm
Standard Foot Shell (With CAP)		FS
Sandal Toe Foot Shell (No CAP)		FTC-4R-1

Black Spectra Sock and
Heel Wedge included in
each Foot Module delivery.

Foot Shell must be
ordered separately.

*Build height based on size 23cm, 26cm, or 29cm Cat 4 foot module with spectra sock, foot shell and 10mm heel height.

Module selection according to activity and weight									
Weight/ Activity level	0-105lbs 0-48kg	106-140lbs 49-64kg	141-175lbs 65-79kg	176-210lbs 80-95kg	211-245lbs 96-111kg	246-280lbs 112-127kg	281-315lbs 128-143kg	316-350lbs 144-159kg	351-365lbs 160-166kg
Low	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9
High	2	3	4	5	6	7	8	9	-

NOTE: Increasing category numbers indicate increased keel response. FOR BILATERAL PATIENTS: Increase by one category for added stability.



DOWNLOAD PRODUCT SHEET



CONTACT US

ORDER




FEET

Now Available in
EVAQ8®



*Highlander MAX is not PDAC Verified.

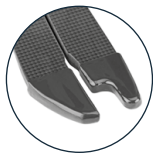

Sandal-toe available for sizes 22-28cm.

Time Tested Stability and Dynamics

The **Highlander®** prosthetic foot provides active patients with long residual limbs an efficient and smooth gait, with excellent energy return. A split heel and keel provides excellent inversion-eversion, allowing users to walk and run with confidence on uneven terrain. The **Highlander MAX** can be ordered to accommodate users up to 500lbs.** As a result, amputees once limited to lower functioning products can now enjoy the benefits provided by a premium energy storage and return device so they can walk farther, faster, longer.

- Water resistant
- Engineered for active amputees with long residual limbs
- Split keel and heel provides excellent inversion-eversion
- Offers ideal blend of stability and action, often preferred by transfemoral amputees
- Can be used for amputees weighing up to 500lbs with the **Highlander MAX** option

TECHNICAL FEATURES

Item N°.	FS3-00			FS3-H5		
Build Height* Highlander	22 - 25cm	3.8" / 97mm	Build Height* Highlander Max	26 - 28cm	4.4" / 111mm	
	26 - 28cm	4.3" / 110mm		29 - 31cm	4.8" / 121mm	
	29 - 31cm	4.6" / 118mm				
Weight (26 Cat 4, module with foot shell and spectra sock) 538g / 1.2lbs			Weight (Cat 10, module with foot shell and spectra sock) 687g / 1.5lbs			
Heel Height			3/8" / 10mm			
Standard Foot Shell (No CAP)			FTC-3M-1			
Standard Foot Shell (With CAP)			FTC-3M-0			
Sandal Toe Foot Shell (No CAP)			FTC-2M-1			

Black Spectra Sock and Heel Bumper Kit included in each Foot Module delivery.

Foot Shell must be ordered separately.

1X998 To avoid spoiling the pyramid during trials



*Build height based on sizes 23, 26, 29, cat. 4 (for Highlander) and cat. 10 (for Highlander MAX), with foot shell, Spectra sock and 10 mm heel height. Sandal-toe available for sizes 22-28cm.

Module selection according to activity and weight													
Weight/ Activity level	100-115lbs 44-52kg	116-130lbs 53-59kg	131-150lbs 60-68kg	151-170lbs 69-77kg	171-195lbs 78-88kg	196-220lbs 89-100kg	221-255lbs 101-116kg	256-285lbs 117-130kg	286-325lbs 131-147kg	326-365lbs 148-166kg	366-420lbs 167-190kg	421-480lbs 191-218kg	481-500lbs 219-227kg
Low	1	1	2	3	4	5	6	7	8	9	10**	11**	12**
Moderate	1	2	3	4	5	6	7	8	9	10**	-	-	-
High	2	3	4	5	6	7	8	9	10**	11**	-	-	-

**Highlander MAX Custom order only




CONTACT US

ORDER




SIERRA®

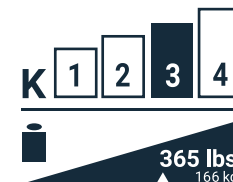

FEET

Now Available in
EVAQ8®

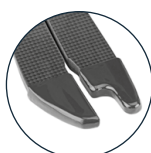

Sandal-toe available
for sizes 22-28cm.

Ultimate Flexibility for a Smoother Gait

The **Sierra®** prosthetic foot incorporates a long carbon fiber spring to provide users a smooth gait over a broad spectrum of activities. The split keel and heel design provides excellent inversion-eversion so that patients can walk with confidence on uneven terrain. A full-length heel enables foot flat early in the gait cycle, providing a compliant foot with a stable base of support. The result is a product that allows amputees to walk comfortably with confidence.

- Water-resistant
- Angle-top design increases the length of the carbon fiber spring, providing a smoother gait and ultimate flexibility

TECHNICAL FEATURES

Item N°.		FS1
Weight (26 Cat 4, module with foot shell and spectra sock)		597g / 1.3lbs
Build Height*	22 - 25cm	5.3" / 135mm
	26 - 28cm	6.2" / 158mm
	29 - 31cm	6.7" / 172mm
Heel Height		3/8" / 10mm
Standard Foot Shell (No CAP)		FTC-3M-1
Standard Foot Shell (With CAP)		FTC-3M-0
Sandal Toe Foot Shell (No CAP)		FTC-2M-1

Black Spectra Sock and
Heel Bumper Kit included in
each Foot Module delivery.

Foot Shell must be
ordered separately.

*Build height based on size 23cm, 26cm, 29cm Cat 4 foot module with foot shell, spectra sock, and 10mm heel height. Sandal-toe available for sizes 22-28cm.

Module selection according to activity and weight										
Weight/ Activity level	100-115lbs 44-52kg	116-130lbs 53-59kg	131-150lbs 60-68kg	151-170lbs 69-77kg	171-195lbs 78-88kg	196-220lbs 89-100kg	221-255lbs 101-116kg	256-285lbs 117-130kg	286-325lbs 131-147kg	326-365lbs 148-166kg
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	-
High	2	3	4	5	6	7	8	9	-	-

1X998 To
avoid spoiling
the pyramid
during trials



CONTACT US


ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com



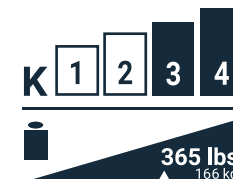

FEET

Now Available in
EVAQ8® and H2O



NOW AVAILABLE!
Sandal Toe
Sizes 22-29cm

Innovation You Can Feel

Thanks to the unique sole plate and virtually indestructible glass composite, the **RUSH RAMPAGE® LP** delivers a smooth transition from heel strike to toe-off. Patients will experience unparalleled comfort, range of motion, and zero dead spots in every step.

- Low-profile build
- Superior ground compliance
- Smooth rollover
- High-energy return
- Cosmesis-friendly
- All-terrain
- Waterproof (fresh or salt water)

TECHNICAL FEATURES

Item N°.		RAMLP
Weight (26 Cat 4, with foot shell and spectra sock)		805g / 1.77lbs
Build Height*	22 - 24cm	4.5" / 114mm
	25 - 27cm	4.7" / 121mm
	28 - 29cm	4.7" / 121mm
Heel Height		3/8" / 10mm
Standard Foot Shell (With CAP)		FS
Sandal Toe Foot Shell (No CAP)		FTC-4R-1

Black Spectra Sock and Heel Wedge included in each Foot Module delivery.

Foot Shell must be ordered separately.

*Build height based on size 23cm, 26cm, or 29cm Cat 4 foot module with spectra sock, foot shell and 10mm heel height.

Module selection according to activity and weight									
Weight/ Activity level	0-105lbs 0-48kg	106-140lbs 49-64kg	141-175lbs 65-79kg	176-210lbs 80-95kg	211-245lbs 96-111kg	246-280lbs 112-127kg	281-315lbs 128-143kg	316-350lbs 144-159kg	351-365lbs 160-166kg
Low	1	2	3	4	5	6	7	8	9
Moderate	2	3	4	5	6	7	8	9	-
High	3	4	5	6	7	8	9	9	-

NOTE: Increasing category numbers indicate increased keel response. FOR BILATERAL PATIENTS: Increase by one category for added stability.

1X998 To avoid spoiling the pyramid during trials



DOWNLOAD PRODUCT SHEET



CONTACT US

ORDER

f i y t #HUMANFIRST



© 2026 by PROTEOR | All rights reserved.

us.proteor.com





ETHOS™ LP

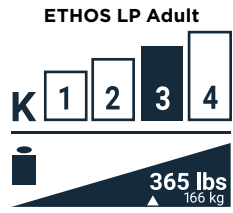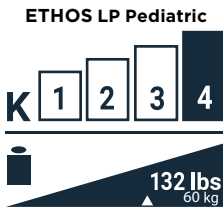

FEET

Available in
EVAQ8®

Low Profile in High Gear

The **ETHOS™ LP** is the ultimate low-profile solution that's integrated with our unique FusionFlex composite technology — a unique layup engineered with a fiberglass core for increased durability and outer layers of carbon fiber for increased energy return. The multiaxial rotation plate between the keel and heel provides independent internal and external rotation to reduce rotational forces on the residual limb in the socket and double the inversion and eversion. This low-profile design offers patients with a longer residual limb higher performance and greater ability to feel better supported in daily activities and achieving goals.

- Durable • High energy return • Smooth rollover • Lightweight • No bolts • Half rocker sole plate • Waterproof

TECHNICAL FEATURES

Item N°.	ETHOS LP for Adults	ETHOS LP for Pediatrics
Weight	485g (26cm Cat 4, module with foot shell)	275g (20cm Cat 2P, module with foot shell)
Build Height*	22 – 25cm = 2.26" / 57mm	19 – 21cm = 1.66" / 42mm
	26 – 28cm = 2.3" / 58mm	
	29 – 30cm = 2.48" / 63mm	
Heel Height	3/8" / 10mm	
Foot Shell	Sandal Toe (No CAP) = FTC-2K-1	Standard Foot Shell (No CAP) = FTC-2M-1

Black Spectra Sock and Heel Wedge included in each Foot Module delivery.
Foot Shell must be ordered separately.

*Build height based on sizes 20cm cat. 2P, 23cm, 26cm, 29cm, cat. 4, with foot shell, spectra sock and 10 mm heel height.

Module selection according to activity and weight - ETHOS LP for Adults										
Weight/ Activity level	100-115lbs 44-52kg	116-130lbs 53-59kg	131-150lbs 60-68kg	151-170lbs 69-77kg	171-195lbs 78-88kg	196-220lbs 89-100kg	221-255lbs 101-116kg	256-285lbs 117-130kg	286-325lbs 131-147kg	326-365lbs 148-166kg
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	-
High	2	3	4	5	6	7	8	9	-	-

Module selection according to activity and weight - ETHOS LP for Pediatrics				
Weight/ Activity level	<58lbs <26kg	59-79lbs 26-36kg	80-100lbs 36-45kg	101-132lbs 46-60kg
High	1P	2P	3P	4P

NOTE: Increasing category numbers indicate increased keel response.



Available in Sandal Toe
only for sizes 22-30cm.

Regular Toe
only for 19-21cm.

DOWNLOAD PRODUCT SHEET




CONTACT US

ORDER



PACIFICA® LP

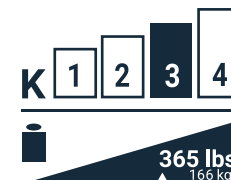


FEET

Because Height Matters

The **Pacifica® LP** prosthetic foot is ideal for users with long residual limbs, due to its low profile design. In contrast with competing low profile foot products, which may feel stiff and uncomfortable, **Pacifica LP** users experience flexibility and energy return allowing them to walk further, faster, and longer. The lightweight and durable design requires less energy expenditure so that amputees are comfortable on their feet for longer periods of time. The pyramid lies within the foot shell, simplifying the cosmetic finishing process. A low build height also provides room for additional componentry, such as a rotator or shock absorber which can provide added comfort.

- Water-resistant
- Perfect for with long residual limbs or those needing more clearance
- LP design requires as little as 2" of clearance
- Lightweight and durable

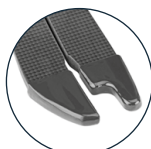

Sandal-toe available for sizes 22-28cm.

TECHNICAL FEATURES

Item N°.		FS4: Pacifica® LP (sizes: 22 - 31cm)
Weight (26 Cat 4, module with foot shell and spectra sock)		537g / 1.1lbs
Build Height*	22 - 25cm	1.9" / 48mm
	26 - 28cm	2.1" / 54mm
	29 - 31cm	2.3" / 58mm
Heel Height		3/8" / 10mm
Standard Foot Shell (No CAP)		FTC-3M-1
Standard Foot Shell (With CAP)		FTC-3M-0
Sandal Toe Foot Shell (No CAP)		FTC-2M-1

Black Spectra Sock and Heel Bumper Kit included in each Foot Module delivery.

Foot Shell must be ordered separately.

*Build height based on size 23cm, 26cm, 29cm Cat 4 foot module with foot shell, spectra sock, and 10mm heel height. Sandal-toe available for sizes 22-28cm.

Module selection according to activity and weight										
Weight/ Activity level	100-115lbs 44-52kg	116-130lbs 53-59kg	131-150lbs 60-68kg	151-170lbs 69-77kg	171-195lbs 78-88kg	196-220lbs 89-100kg	221-255lbs 101-116kg	256-285lbs 117-130kg	286-325lbs 131-147kg	326-365lbs 148-166kg
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	-
High	2	3	4	5	6	7	8	9	-	-

1X998 To avoid spoiling the pyramid during trials



CONTACT US



ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com



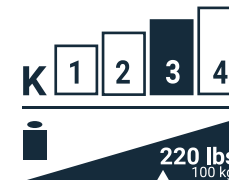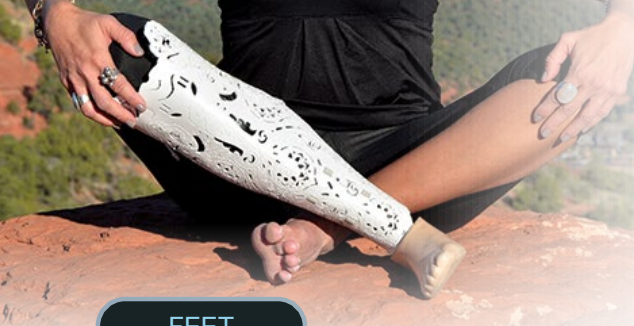

FEET

Comfort and Dynamics

Aided by the shape of the carbon strut, the **DynaStep** ensures energy savings and comfort while providing the patient with dynamics suitable to a K3 activity level; the split keel design allows inversion and eversion while navigating uneven terrain.

- Appropriate for persons with transtibial and transfemoral amputations



PDAC
VERIFIED

L5981



TECHNICAL FEATURES

Item N°.	Side	Size (cm)	Max. Patient Weight
1A101-D...	Right	22 to 25	176lbs / 80kg
1A101-G...	Left	26 to 29	220lbs / 100kg

DynaStep Foot + Multi-axis Ankle Kit

The combination of a **DynaStep** foot and a 1D111 Multiaxis Ankle allows the patient to perform motions of plantar and dorsiflexion, as well as inversion-eversion. Particularly appropriate for users.



PDAC
VERIFIED

L5981 & L5986



TECHNICAL FEATURES

KIT	Includes	Size (cm)
1E201	1A101 + 1D111	22 to 29



Ankle Item N°.

1D111



CONTACT US

ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com





FEET

Comfort and Freedom

The **Gery** foot aims to provide lower level K2 users with the freedom they deserve. It greatly assists your patient's ability to ambulate with its lightweight design and accommodates heel heights ranging from .2" to .8" (5mm to 20mm). The **Gery** is easily cleaned with water and resistant to dirt thanks to its closed cell EVA foam.

- Engineered specifically for the low-level K2 ambulators
- Appropriate for both transtibial or transfemoral amputations
- Not recommended for persons with bilateral amputations
- Lightweight thanks to its EVA structure $\approx$ 30% lighter than a SACH foot!
- 1 Year warranty



TECHNICAL FEATURES

Weight	395g / 0.9lbs
Build Height	3.2" / 82mm
Heel Height (T)	from .2" to .8" (5mm to 20mm)
Maximum patient weight	330lbs / 150kg
Activity level	K2

Item N°.	Side	Size	Color
1A200-G...H	G = Left	22 to 29cm	H = Light
1A200-D...H	D = Right		
1A200-G...BH	G = Left		BH = Dark
1A200-G...BH	D = Right		



DOWNLOAD PRODUCT SHEET



CONTACT US

ORDER

[f](#) [i](#) [v](#) [in](#) #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com





FEET

In the gym, at the park, or on the trail... the EasyRun blade by PROTEOR is the perfect foot to get you moving.

So, you're not an experienced sprinter? Not to worry! Whether you're exploring new activities or looking to take your daily workout to a new level, the **EasyRun** will help you increase strength, mobility, happiness, and amplify your life. With its all-terrain, all-weather, curved fiberglass blade design, and anti-slip Vibram® sole you can take your **EasyRun** workout outside with confidence.

- Rotatable pyramid adapter
- Curved fiber glass blade
- Anti-slip sole.
- 6 vibrant designs



USA



Orange

Black

Water

Tag

Purple

TECHNICAL FEATURES

Weight	2.40lbs / 1088g
Build Height	8.6in / 220mm
Materials	Fiberglass + rubber
Connection	Rotatable Pyramid

Design	Orange	Black	Water	Tag	Purple	USA
Item N°.	1A920-...O	1A920-...N	1A920-...E	1A920-...G	1A920-...V	1A920-...U

... = Module

Green sole kit	1A92013-1
Black sole kit	1A92013-2

Module selection according to activity and weight							
Activity/Weight	88-108lbs 40-49kg	110-130lbs 50-59kg	132-152lbs 60-69kg	154-174lbs 70-79kg	176-194lbs 80-88kg	198-218lbs 90-99kg	220-242lbs 100-110kg
Moderate Impact	-	1	2	3	4	5	6
High Impact	1	2	3	4	5	6	-

1X998 To avoid spoiling the pyramid during trials



DOWNLOAD PRODUCT SHEET



CONTACT US


ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com



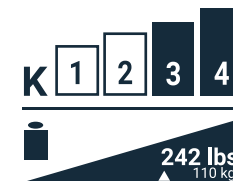

FEET

More versatile than a typical running blade. More dynamic than a sport foot.

Engineered with the highest quality materials for the utmost versatility, durability, and comfort.

ELEVATED BY AIRBUS®

- **Lightweight, ultra-flexible, dynamic carbon fiber blade** is upcycled from higher-quality aerospace-grade material
- **Five levels of stiffness** for diverse types/levels of activity
- **Returns energy** smoothly and comfortably
- **C-shape** allows heel-to-toe movement and rolling motion, like a natural running movement
- **Built to last**, maintains shape over time

ENHANCED BY SALOMON®

- Patented sole specifically designed to **grip like a shoe** is optimal for all terrains
- Enhances **the step through**, providing a **natural stride** feeling
- Maximizes ground contact surface, **improving stability**
- **Interchangeable**, easy to remove/clean/change



Material	Recycled carbon fiber from the aerospace industry
Weight	- 840g (blade + sole + adaptor) for the Audacieuse - 740g for the Intrepide
Versions	- 5 stiffness ranges (XS to XL) for the Audacieuse - 5 stiffness ranges (PXS to PXL) for the Intrepide
Water	Usable in contact with fresh, salt, and chlorinated water.
Outsole	Salomon® Blue Fire outsole (all road) Ref. HS01B
Maximum User Weight	110kg / 242lbs
Minimum Clearance (distance between the bottom of the socket and the ground on the walking prosthesis with shoe)	> 24cm for the Audacieuse > 19cm for the Intrepide

Module selection according to activity and weight											
Weight	77-100lbs 35-45kg	101-121 lbs 46-55 kg	122-148 lbs 56-67 kg	149-172 lbs 68-78 kg	173-199 lbs 79-90 kg	200-242 lbs 91-110 kg					
Activity Level	1/2	3	1/2	3	1/2	3	1/2	3	1/2	3	1/2
Audacieuse: >24cm*	-	-	XS	S	S	M	M	L	L	XL	XL
Intrepide: >19cm*	PXS	PXS	PXS	PS	PS	PM	PM	PL	PL	PXL	PXL

Activity level: 1=Low / 2=Medium / 3=High

FOR BILATERAL PATIENTS: Increase by one category for added stability.



***Minimum clearance:** distance between the ground and the bottom of the socket on the walking prosthesis with shoe.

DOWNLOAD PRODUCT SHEET



CONTACT US



ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com



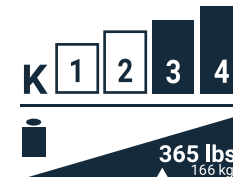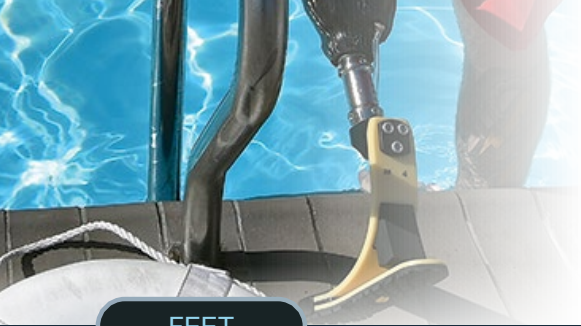

FEET

Now Available in
EVAQ8® and H2O



Incomparable Grip. No Foot Shell Required.

The **RUSH Foot H2O Collection** is the ideal foot for your shower leg or various, non invasive water sport activities. Available in the tried-and-true RUSH Foot **HiPro**, **RAMPAGE**, **RAMPAGE LP** and **ROGUE 2** models and featuring an integrated Vibram® sole and custom molded rubber toe wrap. The **RUSH Foot H2O Collection** of feet offers the same great durability and comfort of all the other RUSH feet, but now adds superior traction on slick surfaces. From the shower to the pool, you won't find a foot that provides this amount of grip and added confidence on slippery surfaces.

- Showering
- Surfing
- Stand Up Paddleboard
- Indoor Yoga
- Gymnastics
- Pool Activities
- Kayaking
- Boarding
- Karate
- Barefoot Activities
- Water proof (fresh and saltwater)

TECHNICAL FEATURES

Item N°.		HiPro: H2H	RAMPAGE LP: H2RAMLP	ROGUE 2: H2R2	RAMPAGE: H2RAM
Weight (26 Cat 4, module only)		765g / 1.7lbs	650g / 1.4lbs	883g / 1.9lbs	695g / 1.5lbs
Build Height*	22 - 24	6" / 152mm	4.1" / 105mm	5.6" / 142mm	6.2" / 159mm
	25 - 27	6.2" / 159mm	4.1" / 105mm	6" / 152mm	6.5" / 165mm
	28 - 29	6.8" / 175mm	4.3" / 111mm	6.5" / 165mm	7" / 178mm

*Build height based on size 23cm, 26cm, or 29cm Cat 4 foot module only.

Black Spectra Sock and
NO heel wedge, each Foot
Module delivery.
No footshells

1X998 To
avoid spoiling
the pyramid
during trials



DOWNLOAD PRODUCT SHEET



CONTACT US

ORDER

f i y t #HUMANFIRST


© 2026 by PROTEOR | All rights reserved.

us.proteor.com





RUSH^{FOOT} EVAQ8 COLLECTION

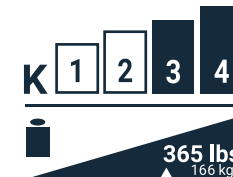

FEET

Now Available in
EVAQ8[®] and H2O



L5987 & L5984
(ROGUE 2)



NOW AVAILABLE!

Sandal Toe for
All EVAQ8 Feet
Sizes 22-29cm

Simplistic, Effective, Efficient

Finally, a vacuum system that seamlessly integrates effective vacuum with the same rugged simplicity of all the RUSH Foot designs. The **RUSH EVAQ8[®] Collection** is the easy solution to providing maximum vacuum efficiency in just a few quick steps. No batteries or complicated mechanical pump required.

- Simple, integrated design
- Stabilizes volume, increases linkage, and heightens proprioception through elevated vacuum suspension
- All-terrain, fully submersable
- Requires minimal maintenance
- The **RUSH EVAQ8 Collection** is available in the **RAMPAGE, HiPro, ROGUE 2, and RAMPAGE LP** models
- All models now available in sandal toe

TECHNICAL FEATURES

Item N°.		RAMPAGE: EVRAM	HiPro: EVQH	RAMPAGE LP: EVGRAMLP	ROGUE 2: EVQR2
Weight (26 Cat 4, with foot shell and spectra sock)		860g / 1.90 lbs	930g / 2.05 lbs	815g / 1.80 lbs	1048g / 2.31 lbs
Build Height*	22 - 24cm	6.6" / 168mm	6.5" / 165mm	4.5" / 114mm	6.1" / 155mm
	25 - 27cm	7.1" / 181mm	6.8" / 175mm	4.7" / 121mm	6.6" / 168mm
	28 - 29cm	7.2" / 184mm	7.1" / 181mm	4.7" / 121mm	6.7" / 171mm
Heel Height		3/8" / 10mm			
See Foot Module Pages for Foot and category					

*Build height based on size 23cm, 26cm, or 29cm Cat 4 foot module with spectra sock, foot shell and 10mm heel height.

Black Spectra Sock and
Heel Wedge included in
each Foot Module delivery.
Foot Shell must be
ordered separately.

REBUILD KIT FOR THE RUSH EVAQ8 COLLECTION = ITEM #EVRB

- Tygon Vacuum Hose - 2 ft.
- Right Angle Hose Adapter
- 90 Degree Hose Barb
- Hose Barb Filter
- Hose Retainer, EVAQ8
- Hose Retainer Clip, EVAQ8
- Straight Hose Barb
- Breather Vent
- Duckbill Valve - Intake
- Duckbill Valve - Exhaust
- EVAQ8 Small O-ring
- EVAQ8 Large O-ring

1X998 To
avoid spoiling
the pyramid
during trials



DOWNLOAD PRODUCT SHEET



CONTACT US

ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com





EVAQ8®

EVAQ8® - The Collection has Expanded!

Meeting the need for innovation, versatility, and pure satisfaction, the popular **EVAQ8** vacuum system collection by PROTEOR is expanding beyond our legacy RUSH feet to carbon fiber, EnduraCore feet, and the newest addition — FusionFlex.

- Optimal design - No additional build height and very low added weight (40g)
- High durability - Rubber vacuum chamber is designed to last - even with continued use over time
- Responds automatically - Pulls vacuum as needed with every step
- The **EVAQ8 Collection** now available in **Shockwave, Agilix, Sierra, Highlander, DynAdapt, Kinterra, and ETHOS LP**



TECHNICAL FEATURES

Item N°.	SHOCKWAVE: F20-V2	AGILIX: F15-V2	SIERRA: FS1-V2	HIGHLANDER: FS3-V2	HIGHLANDER MAX: FS3-V2*	DYNADAPT: F10-V2	KINTERRA: RM3-V2	ETHOS LP: LP3-V3
Weight (26 Cat 4, with foot shell and spectra sock)	1018g 2.24 lbs	645g 1.42 lbs	637g 1.40 lbs	578g 1.27 lbs	687g 1.5 lbs	668g 1.47 lbs	873g 1.92 lbs	485g 1.1 lbs
Build Height*	22-25cm**	6.9" / 176mm	5.5" / 140mm	5.3" / 135mm	3.8" / 95mm	-	6.1" / 155mm	4.56" / 118mm
	26-28cm	7.4" / 189mm	6.2" / 159mm	6.2" / 158mm	4.3" / 110mm	4.44" / 111mm	6.5" / 166mm	4.67" / 119mm
	29-31cm***	7.8" / 200mm	6.8" / 173mm	6.7" / 172mm	4.6" / 118mm	4.6" / 118mm	7.0" / 178mm	4.8" / 122mm
See Foot Module Pages for Foot and category								

*Build height based on size 23cm, 26cm, or 29cm cat. 4 foot module with spectra sock, foot shell, and 10mm heel height. Highlander MAX weight and build height based on size 26cm, 29cm, cat. 10. with spectra sock, foot shell, and 10mm heel height.

**Shockwave = 23-25cm

***ShockWave, Kinterra, and ETHOS LP = 29-30cm

Black Spectra Sock and Heel Bumper kit included in each Foot Module delivery.

Foot Shell must be ordered separately.

REBUILD KIT FOR THE EVAQ8 COLLECTION MODULES ABOVE = ITEM #EV2RB

- Tygon Vacuum Hose – 2 ft.
- Right Angle Hose Adapter
- 90 Degree Hose Barb
- Hose Barb Filter
- Hose Retainer, EVAQ8
- Hose Retainer Clip, EVAQ8
- Straight Hose Barb
- Breather Vent
- Duckbill Valve – Intake
- EV2 Exhaust Valve Spacer

1X998 To avoid spoiling the pyramid during trials



DOWNLOAD PRODUCT SHEET



CONTACT US

ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com



FOOT SHELLS & SPECTRA SOCKS

FEET

PROTEOR offers a variety of cosmetic foot shells and Spectra Socks.

Please use the following guide to determine which foot shells and socks are available for each foot product. Not all Foot Modules are available in Sandal Toe; ensure the corresponding Sandal Toe foot module is also ordered (CAPs are not available for sandal toe products).

Carbon Fiber / EnduraCore / FusionFlex Modules		___-___-__	___	___	___-__	___
Product	Foot Shell Type	Product ID	Size	Shell Color L: Light / M: Medium D: Dark / B: Jet Black	Toe Configuration	Side L: Left R: Right
PROTEOR Agilix* DynAdapt* Highlander* Highlander MAX* Pacifica LP Sierra*	Regular Toe, with Cap	FTC-3M-0	22-31	L / M / D	4-R	L / R
	Regular Toe, No Cap	FTC-3M-1	22-31	L / M / D / B	4-R	L / R
	Sandal Toe, No Cap	FTC-2M-1	22-28	L / M / D / B	4-S	L / R
ETHOS LP PEDS	Regular Toe, No Cap	FTC-2M-1	19-21	L / M / D / B	4-R	L / R
ETHOS LP*	Sandal Toe, No Cap	FTC-2K-1	22-30	L / M / D / B	4-S	L / R
ShockWave*	Regular Toe, No Cap	FTC-3F-1	23-30	L / M / D / B	4-R	L / R
	Sandal Toe, No Cap	FTC-2F-1	23-28	L / M / D / B	4-S	L / R
PROTEOR Kinnex 2.0	Regular Toe, No Cap	FTC-3K-1	24-30	L / M / D / B	4-R	L / R
	Sandal Toe, No Cap	FTC-2K-1	24-28	L / M / D / B	4-S	L / R
Kinterra*	Sandal Toe, No Cap	FTC-2K-1	22-30	L / M / D / B	4-S	L / R

*Includes EVAQ8 foot modules

**Foot shell compatibility for the older Kinterra 2.0 module: FTC-3M and FTC-2M

FOOT SHELL PART NUMBER EXAMPLE: To order an Agilix Foot Shell in a size 26cm right-side sandal toe, medium skin tone, the part number would be: FTC-2M-126M4-SR



CONTACT US

ORDER

#HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com



FOOT SHELLS & SPECTRA SOCKS

FEET

PROTEOR offers a variety of cosmetic foot shells and Spectra Socks.

Please use the following guide to determine which foot shells and socks are available for each foot product. Not all Foot Modules are available in Sandal Toe; ensure the corresponding Sandal Toe foot module is also ordered (CAPs are not available for sandal toe products).

RUSH — Sandal Toe Foot Shell		— — — — —	— —	—	— — —	—
Product	Foot Shell Type	Product ID	Size	Shell Color L: Light / M: Medium D: Dark / B: Jet Black	Toe Configuration	Side L: Left R: Right
RUSH ROGUE 2* RUSH HiPro* RUSH RAMPAGE* RUSH RAMPAGE LP*	Sandal Toe, No Cap	FTC-4R-1	22-29	L / M / D / B	4-S	L / R

RUSH — Regular Toe Foot Shell		— — — — —	— — — — —	— — — — —	—
Product	Foot Shell Type	Product ID	Size	Side L: Left R: Right	Shell Color 1: Light 3: Dark
RUSH ROGUE 2* RUSH HiPro* RUSH RAMPAGE* RUSH RAMPAGE LP*	Regular Toe, with Cap	FS	22-29	L / R	1 / 3
RUSH Chopart	Regular Toe, with Cap	FS	22-30	L / R	1 / 3
RUSH Chopart PEDS	Regular Toe, with Cap	FSK	19-21	L / R	1 / 3

Spectra Socks



Description	Item N°.
Black Spectra Sock (1 pk)	
22-25cm foot module = 30cm/12"	SO-NPS-200XX-00 (example: size 26cm foot module = SO-NPS-20035-00)
26-28cm foot module = 35cm/14"	
29-30cm foot module = 40cm/16"	
Kinnex White Custom Spectra Sock (1 pk)	
Sizes: 24 - 30cm	SO-F14-000XX-00 (example: size 26cm foot module = SO-F14-00026-00)

*Includes EVAQ8 foot modules

FOOT SHELL PART NUMBER EXAMPLE: To order a RUSH ROGUE 2 Foot Shell in a size 26cm right-side sandal toe, medium skin tone, the part number would be: FTC-4R-126M4-SR



CONTACT US

ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com





Digital


Great outcomes are built in layers: a scan, a design, a print, a final shape. Each step turns a little more data into something a patient can actually wear.

From scan to fabrication, the PROTEOR Digital Suite connects scanning, modeling, 3D printing, and milling into one workflow. The result is a device built with accuracy, speed, and craftsmanship at every stage.



Our full packages

have been specially designed to support you from A to Z in the design and fabrication of your devices!

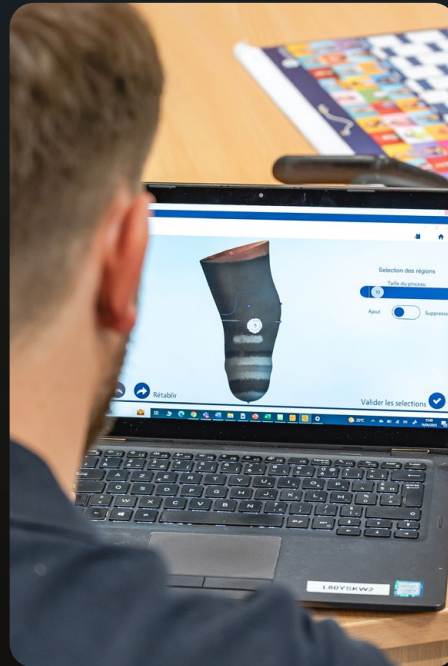


[CONTACT US](#)

[ORDER](#)

Discover the different ways of starting the working flow with your patient:

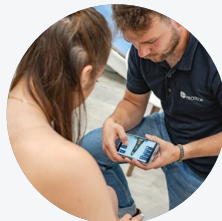

The patient cannot be moved	The patient can sit, stand or be moved	
Model function in the « Proteor Scan » app	2 methods depending on devices needs available in the « Proteor Scan » app	
MODEL 	STRUCTURE SENSOR 	TRUEDEPTH 
Select the Model function and enter your measurements directly	Scan with the Structure Sensor scanner	Scan with Apple's Truedepth technology (Face ID)
Emergency brace, neck brace, AFOs/KAFOs, seating device, trans-femoral prosthesis	Scoliosis brace, immobilizing brace, seating device, standing device, orthosis	Lower and upper limb prosthesis
 More respectful	 Fast	 Color scanning

- Free smartphone application
- Available on iOS (iPad and iPhone)
- Cloud backup
- Works offline (last 30 scans stored locally)
- Compatible with Structure Sensor: Original, Mark II, Pro, 3
- Ordering information: ref 8A09 - Scanner structure Sensor 3 (without cable)
- Languages available: chinese, czech, english, french, german, italian, japanese, portuguese, russian, spanish
- Structure Sensor warranty: 1 year




CONTACT US

ORDER

Advance your practice with an intuitive & complete software. 2 solutions to meet your needs:

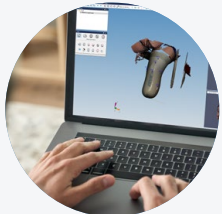


FREESTYLE DESIGN	STEP BY STEP DESIGN
	
	
Wide toolbox to modelize and design the medical device	Structured workflow
Creative freedom Complex designs without restrictions	Guided steps Easy to assimilate
Computer-based modeling	iPad / iPhone or Computer-based modeling

- Available on PC Windows 10-11
- Interconnected with « Proteor Scan », « Proteor Carve » and « Proteor Print »
- Open file formats: stl, obj
- Languages available: chinese, czech, english, french, german, italian, japanese, portuguese, russian, spanish
- Cloud backup



VOLUME-BASED PRICING

Payment per export file, whatever production solution

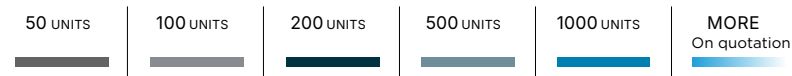

* 1 unit= 1 file export to mill or print.

CONTACT US

ORDER

Introduce the latest 3D printing technologies to your production tools for orthopedic devices.



OPEN PRODUCTION

3D PRINTERS	ICARUS PRO	ICARUS LITE
Item number	PRO3DP.F	LITE3DP.F
Custom-made orthopedic devices	Test socket Definitive socket Trunk Brace Others ...	Test socket Definitive socket Flexible inner Others...
Material	CPX-KM, PETG, TPU, PP, EVA+, EVA Form	
Software compatibility	Proteor Shape and Others	
Dimensions in m (Width x Depth x Height)	1,346 m wide x 0,813 m deep x 1,702 m high	0,660 x 0,760 x 1,240
Maximum print size in mm (Width x Depth x Height)	1,346 x 0,813 x 1,702 and 700 x 500 x 1000	305 x 305 x 610
Weight in kg	225	82
Filament diameter in mm	2,85	
Average printing time	Test Socket: between 1h30 to 2h00 Definitive Socket: between 4h00 to 8h00	

MATERIALS

Item number	Printers model	Material	Weight kg	Ø wire mm	Color
CPX.KM.F.2.85.3.2kg	ICARUS / ICARUS LITE	CPX (Polypropylene + carbone)	3,2	2,85	black
PETG.NAT.F.2.85.4.5kg	ICARUS / ICARUS LITE	PETG	4,5	2,85	natural white
CPX.PP.NAT.F.2.85.3.5kg	ICARUS / ICARUS LITE	Polypropylene	3,5	2,85	natural white
EVAF-285-3.5KG-BLK	ICARUS / ICARUS LITE	EVA Form	3,5	2,85	black
EVAP-285-3.5KG-NAT	ICARUS / ICARUS LITE	EVA +	3,5	2,85	natural white

Note: EVA+ and EVA Form come in alternate colors as well.

CONTACT US

ORDER

Carving solutions meeting your orthopedic needs

Three available models (Classic/Complete/Premium) to meet different device needs, allowing you to rapidly mill the positive rectified shape of a patient in a polyurethane foam block with excellent surface quality.



Time-tested, all-day milling

- All three models have been designed with industrial conception to facilitate full-day milling.
- Same carver models used daily, over many years by PROTEOR Central Fabrication.



Extensive expertise

- PROTEOR's internal experience as fabricators ensures the machine is set with optimized parameters to provide the best milling result.
- Optimal surface quality, limiting manual operations after milling and optimizing time.



Integration ease

- PROTEOR carvers can be easily integrated into your workshop.
- We meet with you first to discuss your individual needs before installation.
- Compatible with many software programs (not just PROTEOR Shape).
- Our support team is readily available to facilitate a seamless transition.

Models	PROTEOR Carve Complete	PROTEOR Carve Classic	PROTEOR Carve Premium
Nb of axis	Yes	No	No
Half-axis	Yes	No	No
Undercut	Yes	No	No
Maximum carved volume (mm)	H: 1000 Ø : 800	H: 1000 Ø : 800	H: 1000 Ø : 800 On the table : L : 1000 W : ± 400 T : 380
Example of orthopedic devices	Spinal orthosis AK / BK socket A.F.O. / K.A.F.O. Cranial remolding orthosis	Spinal orthosis AK / BK socket	Spinal orthosis AK / BK socket Foam bed Verticalizer Seat brace
Machine dimensions	Width : 1860mm Depth : 1560mm Height : 2005mm	Width : 1860mm Depth : 1560mm Height : 2005mm	Width : 2080mm Depth : 2350mm Height : 2000mm

CONTACT US

ORDER

Connectors

Strength is often found in what holds things together. In alignment. In continuity. In the assurance that each part supports the next.

Connectors create integration and stability across the system. Energy transfers smoothly from one element to the next, maintaining alignment through every shift in weight and direction. Separate components function as one continuous experience.



CONNECTORS & COMPONENTS

TABLE OF CONTENTS



KITS



PYLON



ADAPTERS



LAMINATION ADAPTERS



CONNECTORS



4-HOLE

CONTACT US

ORDER

PYLONS

COMPONENTS

Kits


	Item N°.	Description	Max patient weight	Material	Weight
	TK200	Aluminum 10" modular 3 piece set. Includes 30mm Aluminum pylon with receiver, tube clamp with water resistant fasteners and 4-hole pyramid with through hole	220lbs / 100kg	Aluminum	-
	TK300	Titanium 10" modular 3 piece set. Includes 30mm Aluminum pylon with receiver, tube clamp with water resistant fasteners and 4-hole pyramid with through hole	365lbs / 165kg	Titanium	395g
	KIT-00-20200-00	Heavy Duty Adapter Kit for Highlander MAX. Includes (2) 34mm tube clamps, (1) 4-ear connector for 34mm system and (1) carbon pylon/ tube	500lbs / 227kg	Titanium and Graphite	744g

[CONTACT US](#)
[ORDER](#)

PYLONS

COMPONENTS

Pylons & Pylon Clamps

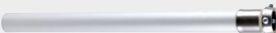	Item N°	Description	Max patient weight	Material	Weight	Total Height	Build Height
	1D52	Pylon Ø 30mm with Pyramid Receiver	220lbs / 100kg	Aluminum	111g	130mm	63mm
	ACC-00-20203-00	Heavy Duty Ø 34mm - 305mm length Pylon/Tube for Highlander MAX	500lbs / 227kg	Graphite	269g	305mm	-
	1D52-HD	Pylon Ø 34mm with Pyramid Receiver	275lbs / 125kg	Aluminum	121g	130mm	67mm
	1D45	Pylon adapter Ø 30mm with Aluminum Pyramid Receiver	220lbs / 100kg	Aluminum	290g	440mm	440mm
	1D46	Pylon adapter Ø 30mm with Titanium Pyramid Receiver	220lbs / 100kg	Aluminum + Titanium	280g	440mm	440mm
	1G0120#kit	Aluminum pylon Ø 30mm - 200mm length and Pylon clamp	220lbs / 100kg	Aluminum	195g	217mm	-
	1G01-HD	Aluminum pylon Ø 34mm - 420mm length	330lbs / 150kg	Aluminum	235g	420mm	-
	1G07	Aluminum pylon Ø 22mm - 300mm length used for pediatric prosthesis	99lbs / 45kg	Aluminum	105g	300mm	-

Peds

CONTACT US

ORDER

ADAPTERS

COMPONENTS

Tube Clamp Adapters

	Item N°	Description	Max patient weight	Material	Weight	Total Height	Build Height
Peds		1D39 Pediatric Tube Clamp Adapter Used with the Ø 22mm tube for child prosthesis	99lbs / 45kg	Aluminum	45g	38mm	12mm
		ACC-00-20201-00 Tube Clamp 34mm used with Highlander MAX	500lbs / 227kg	Titanium			
		1D40 Tube Clamp Adapter Female 4-screw adapter for tightening on Ø 30mm tube	220lbs / 100kg	Titanium	90g	50mm	16mm
		1D41 Tube Clamp Adapter Female 4-screw adapter for tightening on Ø 30mm tube	220lbs / 100kg	Aluminum	100g	47mm	13mm
		1D41-HD Tube Clamp Adapter for Ø 34mm Tube Female 4-screw adapter for tightening on Ø 34mm tube	330lbs / 150kg	Aluminum	120g	55mm	17mm
		1K182 Threaded Pyramid insert M36 x 1.5mm	275lbs / 125kg	Stainless steel	65g	28mm	6mm
		1K183 Threaded Receiver insert M36 x 1.5mm Use with 1K179 or 1S145	330lbs / 150kg	Stainless steel	75g	23mm	14mm
		1K184 Threaded Pyramid For screwing on 1P200-KD, 1P320-KD, 1P321-KD, 1P360-KD, 1P130-KD, 1P180-KD	275lbs / 125kg	Stainless steel	94g	30mm	15mm

[CONTACT US](#)
[ORDER](#)

CONNECTORS

COMPONENTS

Double Adapters

	Description	Item N°.	Max patient weight	Material	Weight	Total Height	Build Height
	Double Pyramid Receiver.	1K207-032-HD	330lbs / 150kg	Aluminum	100g	32mm	32mm
		1K207-045-HD			130g	45mm	45mm
		1K207-060-HD			145g	60mm	60mm
		1K207-075-HD			160g	75mm	75mm
	Double Pyramid Adapter	1K190		Stainless steel	120g	36mm	7mm

[CONTACT US](#)
[ORDER](#)

CONNECTORS

COMPONENTS

Connectors for PROTEOR Knees

These connectors are necessary for the assembly of PROTEOR knees: 1M01 - 1M01-HD - 1M10 - 1M10-HD - 1M102 - 1M102-HD - 1M102V - 1M102V-HD - 1M112 - 1M113

This specific connection system allows you to select the most appropriate linking piece between the knee and the socket. The simplified shape of the upper part of the knee can receive the proximal linking piece using a M8 screw. The knees listed above must be ordered together with one of the following connectors. All these connectors provide a $\pm 15^\circ$ socket rotation.


	Item N°	Description	Max patient weight	Material	Weight	Connection	Total Height	Build Height	Knee
Peds	 1K160	Pyramid Connector provides the male pyramid connection	220lbs / 100kg	Aluminum + Stainless steel	45g	M8	30mm	7mm	1M01 1M10 1M102 1M102V
	 1K160-HD	Pyramid Connector for 125kg knee provides the male pyramid connection	275lbs / 125kg	Aluminum + Stainless steel	90g			8mm	1M01-HD 1M10-HD 1M102-HD 1M102V-HD
	 1K163	Pyramid Connector provides the male pyramid connection	220lbs / 100kg	Stainless steel	100g				1M01 1M10 1M102 1M102V
	 1K164	Pediatric Pyramid Connector provides the male pyramid connection to the 1M112 knee	99lbs / 45kg	Stainless steel	120g		21mm	9mm	1M112 1M113
	 1K20	This Ø 30mm Tube Clamp provides a long and lightweight link between the knee and a short residual limb socket. To be assembled using Adapter 1K14	220lbs / 100kg	Aluminum + Steel	210g	Tube Ø 30mm M8	218g	210mm	
	 1K30	Aluminum Linking plate is convenient for assembly to the distal end of a transfemoral socket with European 4-hole pattern	220lbs / 100kg	Aluminum	80g	M8	21mm	13mm	

[CONTACT US](#)
[ORDER](#)

LAMINATION ADAPTER

COMPONENTS

Lamination Adapters

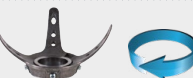	Item N°	Description	Max patient weight	Material	Weight	Total Height	Build Height
	1K40	Used for knee disarticulation or long residual limb M8 connection	275lbs / 125kg	Steel	120g	120mm	6.5mm
	ACC-00-20202-00	4-Ear Connector for 34 mm System used with Highlander MAX	500lbs / 227kg	Titanium	141g	39mm	-
	1K173	3-prong anchor with pyramid	275lbs / 125kg	Stainless steel	150g	65mm	12mm
	1K176	3-prong anchor with rotating receiver	275lbs / 125kg	Aluminum + Stainless steel	205g	68mm	25mm
	1K177	3-prong anchor with rotating pyramid	220lbs / 100kg	Stainless steel	225g	76mm	17mm
	1K179	Transfemoral Threaded base 3-Prong Lamination Adapter for knee disarticulation, delivered with a spacer. To use with M36 Threaded knees or Threaded inserts	330lbs / 150kg	Stainless steel	125g	45mm	12mm
	1K202	4-ear Lamination Anchor	275lbs / 125kg	Titanium	55g	32mm	10mm
	1K203	4-ear Lamination Anchor HD Stainless steel, heavy duty	330lbs / 150kg	Stainless steel	85g	32mm	10mm
	1K204	4-ear Offset Lamination Anchor HD Transfemoral adapter with 9mm axial offset of the pyramid	330lbs / 150kg	Stainless steel	113g	32mm	11mm


[CONTACT US](#)
[ORDER](#)

4-HOLE ADAPTER

COMPONENTS

4-Hole Adapters

		Item N°.	Description	Max patient weight	Material	Weight	Total Height	Build Height
Peds		1K206	Pyramid Adapter for pediatric sockets	99lbs / 45kg	Aluminum + Stainless steel	45g	18mm	6mm
		1K165	Pyramid Receiver for pediatric sockets	99lbs / 45kg	Aluminum	35g	13mm	13mm
		1K172	Stainless Steel Pyramid Receiver for Socket adapter features European 4-hole pattern and a pyramid receiver allowing for tilt adjustment of the socket	275lbs / 125kg	Stainless steel	100g	16mm	16mm
	 	1K66	4-hole Receiver Socket adapter features European 4-hole pattern and a pyramid receiver allowing for tilt adjustment of the socket and rotation adjustment	275lbs / 125kg	Titanium	93g	26mm	26mm
		1K205	4-hole Pyramid Socket adapter features European 4-hole pattern and a pyramid adapter allowing for tilt adjustment of the socket	275lbs / 125kg	Stainless steel	80g	22mm	22mm
	 	1K60	4-hole Pyramid Socket adapter features European 4-hole pattern and a pyramid adapter allowing for tilt adjustment of the socket and rotation adjustment	275lbs / 125kg	Titanium	106g	23mm	23mm
		1K209	4-hole Pyramid Socket Adapter features European 4-hole pattern and a pyramid adapter allowing for tilt adjustment of the socket	275lbs / 125kg	Titanium	46g	22mm	8mm
		1K14	Ø 30mm Tubular Link providing connection between a 30mm Ø tube and a socket with a 4-hole plate	220lbs / 100kg	Aluminum	60g	32mm	32mm

[CONTACT US](#)
[ORDER](#)

4-HOLE ADAPTER

COMPONENTS

4-Hole Lamination Adapters

	Use	Material	Weight	Max. Patient Weight	Item N°.
	Threaded Four Hole Lamination Adapter With One-Way Valve A threaded 4-hole adapter plate with one-way valve. Kit includes #10 plate, TP5 valve, and four attachment screws	Aluminum	103.6g	300lbs 136kg	FHLA-10K
	Threaded Four Hole Lamination Adapter With One-Way Valve A threaded 4-hole adapter plate with one-way valve. Kit includes #10 plate, TP3 valve, and four attachment screws	Aluminum	103.6g	300lbs 136kg	FHLA-10KT-P3SB
	Threaded Four Hole Lamination Adapter With One-Way Valve A threaded 4-hole adapter plate with one-way valve. Kit includes #16 plate; 1 - TP7 & 1 - TP10 valve	Aluminum	93.2g	300lbs 136kg	FHLA-16K
	Threaded Four Hole Lamination Adapter A threaded 4-hole adapter plate. Compatible with TP9, TP5 and TP3 valves	Aluminum	98.9g	300lbs 136kg	FHLA-10

[CONTACT US](#)
[ORDER](#)

4-HOLE ADAPTER

COMPONENTS

4-Hole Lamination Adapters

	Use	Material	Weight	Max. Patient Weight	Item N°.
	Threaded Four Hole Lamination Adapter Standard 65mm hole pattern. Dimensions: 72 x 11mm	Aluminum	60.96g	300lbs 136kg	FHLA-11
	Threaded Four Hole Lamination Adapter A low-profile, threaded 4-hole adapter plate. Dimensions: 72 x 9.5mm Compatible with TP9, TP5 and TP3 valves	Aluminum	53g	300lbs 136kg	FHLA-13
	Threaded Four Hole Lamination Adapter Standard 65mm hole pattern. Dimensions: 69 x 9.5mm	Aluminum	56.7g	300lbs 136kg	FHLA-3
	Threaded Four Hole Lamination Adapter M6 threaded blind holes. Plate has no through holes. Dimensions: 72 x 15mm	Nylon	56.7g	300lbs 136kg	FHLA-23
	Threaded Four Hole Lamination Adapter M6 threaded blind holes. Plate has no through holes. Dimensions: 72 x 15mm	Nylon	60g	300lbs 136kg	FHLA-24

CONTACT US

ORDER



Socket Technology

Confidence starts with comfort. When something fits well, it becomes part of us. It allows focus to shift away from mechanics and back to living.

Socket technology distributes pressure thoughtfully and reinforces stability for sustained wear. The connection between body and system feels secure and dependable, allowing focus to return to daily life.

SOCKET TECHNOLOGY

TABLE OF CONTENTS

Keasy



PREFAB CONE

STEPLINE PLUS



SEALING SLEEVES



SEALING RING



ADJUSTABLE TEMPORARY SOCKET



STEP-IN

SILITOUCH



THE COMPRESSION ADVANTAGE



LOCKS AND VALVES

- Valve for thermoformable transfemoral socket
- Valve for Aluminum transfemoral socket
- Accessories for transfemoral valve
- Depressurizing valves
- Locks for distal attachment
- Tools and accessories for lock



COSMETICS

- Transfemoral 2-part cosmetic cover
- Polyethylene & Polyurethane one-part cosmetic cover
- Plastazote one-part cosmetic cover

CONTACT US

ORDER

LINERS


SOCKETS



KEASY has established itself over the years as the unique prefabricated cone allowing to manufacturing rapidly quality finished thermoformed liners.

A unique prefabricated cone, the **KEASY** has established itself due to the ability to rapidly create quality, finished thermoformed liners.

No welding, no gluing, and seamless, the **KEASY** minimizes allergy risks while providing comfort. Fabrication and finishing time is shortened with easy thermoforming, cutting and sanding abilities. Changes in limb volume can be accommodated by gluing extra material to the outer surface. Closed cell structure keeps perspiration and cleaning products from penetrating the material for easy cleaning.

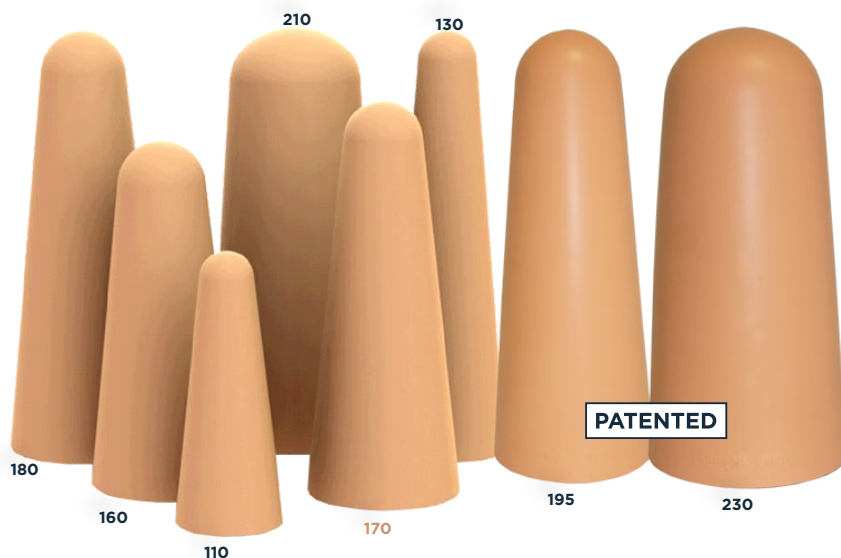

Item N°	Use	Length	Ø Proximal	Profile	Distal Thickness	Packaging
Keasy Comfort						
1S402-110#12	BK small Upper limb	305mm	110mm	5mm	8mm	x 12
1S400-130#24	BK, Long & thin	500mm	130mm	6,5mm	10mm	x 24
1S400-160#24	BK Short & medium	400mm	160mm	6,5mm	10mm	x 24
1S400-180#12	BK Short & medium	500mm	180mm	6,5mm	15mm	x 12
1S400-195#9	BK XL	510mm	195mm	6,5mm	10mm	x 9
1S400-210#6	AK	500mm	210mm	6,5mm	10mm	x 6
1S400-230#6	AK XL	510mm	230mm	6,5mm	10mm	x 6
Keasy interface						
1S401-170#12	BK Short & medium	450mm	170mm	3mm	3mm	x 12



CONTACT US

ORDER



© 2026 by PROTEOR | All rights reserved.

us.proteor.com



LINERS

SOCKETS



HOW TO DETERMINE THE SIZE

Metric measurement of the residual limb circumference, at 4cm from the distal end.

SiLiTOUCH - The Compression Advantage

SiLiTouch is a premium transtibial liner system engineered around consistent compression for a healthier limb. SiLiTouch Classic and SiLiTouch Aloe are built for specific patient needs, with the durability and sleek design clinicians and patients trust. All SiLiTouch liners integrate a 10cm matrix for comfort and suspension.

SiLiTOUCH Classic - Compression First.

SiLiTouch Classic delivers steady, even compression that supports circulation and helps manage limb volume throughout the day, so your patients stay stable, protected, and ready to move in all impact levels. Finished in a sleek, sophisticated black Exterior, it looks as sharp as it performs.

SiLiTOUCH ALOE - Compression, Softened.

SiLiTouch Aloe pairs consistent compression with a softer silicone infused with organic Aloe Vera oil, long trusted in wound care, to protect sensitive or fragile skin. SiLiTouch Aloe is advised for low and moderate impact levels.


Item No. (add size)	Model	Distal thickness (mm)	Thickness (mm)	Sizes to choose from (cm)	Distal attachment
1S712-3-xxx	SiLiTouch - LINER TT Classic locking 3mm - xx	17.5 to 18.5	3	16 - 18 - 20 - 22 - 23.5 - 25 - 26.5 - 28 - 30 - 32 - 34 - 36 - 40 - 45	locking
1S712-6-xxx	SiLiTouch - LINER TT Classic locking 6mm - xx	19 to 20	6	20 - 22 - 23.5 - 25 - 26.5 - 28 - 30 - 32	locking
1S722-3-xxx	SiLiTouch - LINER TT Classic cushion 3mm - xx	14.5 to 16	3	16 - 18 - 20 - 22 - 23.5 - 25 - 26.5 - 28 - 30 - 32 - 34 - 36 - 40	cushion
1S722-6-xxx	SiLiTouch - LINER TT Classic cushion 6mm - xx	16 to 17.5	6	16 - 18 - 20 - 22 - 23.5 - 25 - 26.5 - 28 - 30 - 32 - 34 - 36	cushion
1S711-3-xxx	SiLiTouch - LINER TT Aloe locking 3mm - xx	17.5 to 18.5	3	16 - 18 - 20 - 22 - 23.5 - 25 - 26.5 - 28 - 30 - 32 - 34 - 36 - 40 - 45	locking
1S711-6-xxx	SiLiTouch - LINER TT Aloe locking 6mm - xx	19 to 20	6	16 - 18 - 20 - 22 - 23.5 - 25 - 26.5 - 28 - 30 - 32 - 34 - 36 - 40 - 45	locking
1S721-3-xxx	SiLiTouch - LINER TT Aloe cushion 3mm - xx	14.5 to 16	3	16 - 18 - 20 - 22 - 23.5 - 25 - 26.5 - 28 - 30 - 32 - 34 - 36 - 40	cushion
1S721-6-xxx	SiLiTouch - LINER TT Aloe cushion 6mm - xx	16 to 17.5	6	16 - 18 - 20 - 22 - 23.5 - 25 - 26.5 - 28 - 30 - 32 - 34 - 36	cushion

[CONTACT US](#)
[ORDER](#)

SEALING SLEEVES

SOCKETS

STEPLINE PLUS

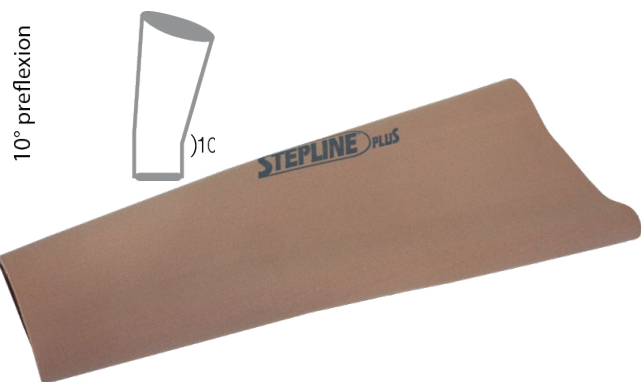

STEPLINE PLUS Sealing Sleeve

The **STEPLINE PLUS** sealing sleeve is made from our copolymer gel and features a durable fabric outer-covering that offers great flexibility. The **STEPLINE PLUS** is available in two colors.



Item N°.	Material	Length	Profile	Distal Circumference	Proximal Circumference	Color
1S500-M	Copolymer Gel	350mm	3mm uniform	21 to 30cm	29 to 43cm	Beige
1S500-L				28 to 48cm	36 to 51cm	
1S500-LP				31 to 46cm	39 to 55cm	
1S500-XL				34 to 50cm	42 to 60cm	
1S500-XXL				36 to 55cm	44 to 70cm	
1S510-M				21 to 30cm	29 to 43cm	Grey
1S510-L				28 to 48cm	36 to 51cm	
1S510-LP				31 to 46cm	39 to 55cm	
1S510-XL				34 to 50cm	42 to 60cm	
1S519-XXL				36 to 55cm	44 to 70cm	

STEPLINE PLUS PREFLEX Sealing Sleeve

Sealing sleeve with 10° pre-flexion. Designed to minimize discomfort during flexion.

Item N°.	Material	Length	Profile	Distal Circumference	Proximal Circumference	Color
1S520-S	Copolymer Gel	340mm	3mm uniform	28 to 48cm	36 to 51cm	Beige
1S520-M				31 to 46cm	39 to 55cm	
1S520-L				34 to 50cm	42 to 60cm	
1S520-XL				36 to 55cm	44 to 70cm	



CONTACT US

ORDER

f i y t #HUMANFIRST

© 2026 by PROTEOR | All rights reserved.

us.proteor.com



SEALING RINGS

SOCKETS

Sealing Ring

This silicone sealing ring has to be bonded onto the liner of lower or upper limb amputees.

By creating vacuum between the socket and the liner, it provides the prosthesis suspension through depressurization for patients without locking device. Available in 4 sizes (from 2 to 5).



TECHNICAL FEATURES

Item N°.	Size	Ø Distal	Adhesive included
1S950-02	2	16 to 22 cm	yes
1S950-03	3	20 to 38 cm	yes
1S950-04	4	30 to 48 cm	yes
1S950-05	5	40 to 55 cm	yes
1S950-12	2	16 to 22 cm	no
1S950-13	3	20 to 38 cm	no
1S950-14	4	30 to 48 cm	no
1S950-15	5	40 to 55 cm	no

ADDITIONAL ITEMS

Item N°.	Description
XC950	Glue for sealing ring
1S950-90	Spray donning Aid for sealing ring

[CONTACT US](#)[ORDER](#)

SOCKETS

SOCKETS

Adjustable Temporary Socket

This socket with adjustable closing allows an easy fitting for patients with recent transfemoral amputation.

- A wide range of adult and child sizes
- A suspension system is necessary (suspenders)
- We strongly recommended to protect the residual limb with cotton tubular jersey



TECHNICAL FEATURES

Item N°.	Side	Item N°. line	Circumference of residual limb base	Height Perineum / Socket End	Lower Connection	Use
1N06-G... or 1N06-D...	G = Left D = Right	3438	34 to 38cm	26.5cm	1N106	Adult
		3943	39 to 43cm			
		4448	44 to 48cm			
		4953	49 to 53cm			
		5458	54 to 58cm			
1N07-G... or 1N07-D...	G = Left D = Right	5963	59 to 63cm	21.5cm	1N107	Child
		3033	30 to 33cm			
		3438	34 to 38cm			

ORDER ONE OF THE FOLLOWING TO COMPLETE YOUR ORDER

	Item N°	Description
	1N106	Linking Device for Temporary Adult Prosthesis
	1N107	Linking Device for Temporary Child Prosthesis
	TJ001/2/3/4/5	Cotton Tubular Jersey
	1X170	Suspenders for Temporary Prosthesis




CONTACT US

ORDER

LINER AID

SOCKETS

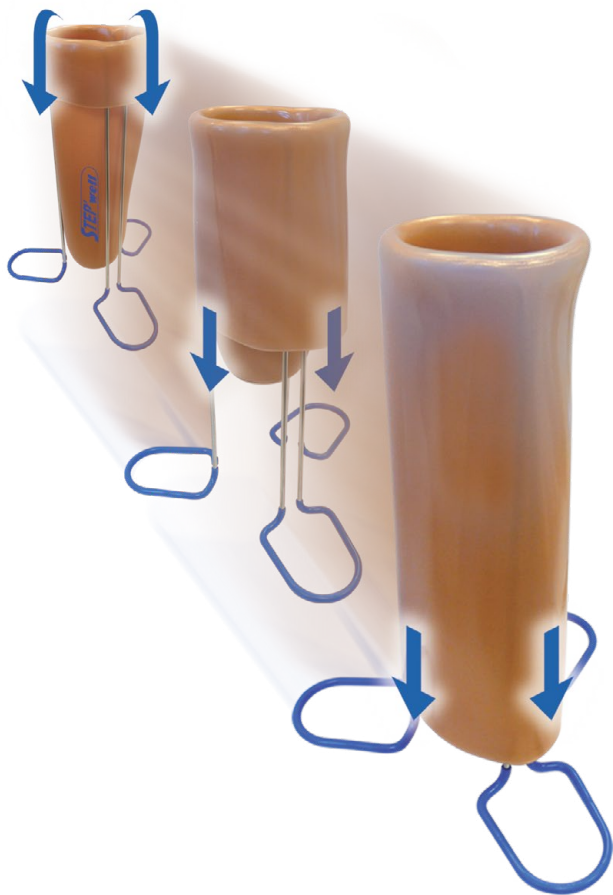

Step-In

Step-in is a multi-purpose support designed to easily reflect and assist in donning liners of any kind (copolymer gel, silicone gel, silicone, etc).

- Facilitates donning, drying, and cleaning of the liner
- Height: 39cm

Item N°.	Packaging
1S900	1
1S900#5	Set of 5

[CONTACT US](#)[ORDER](#)

LOCKS & VALVES

SOCKETS

Expulsion Valves

	Use	Description	Item N°.
	Expulsion Valves for Transtibial Socket Allows a firm suspension of the prosthesis onto the residual limb, due to automatic air expulsion Easy to use for the patients	Valve for BK Socket	1S162
	Expulsion Valves for Transfemoral Socket Allows a firm suspension of the prosthesis onto the residual limb, thanks to manual air expulsion or automatic one. The cap has a quick opening system (1/4 turn), allowing to use a putting on aid for the prosthesis. Easy to use for the patient. A specific tool (1S164) can be used to set up the valve	Valve for AK Socket	1S163
	Valve Tool Allows to set up, remove and configure (manually or automatically) transfemoral valve 1S163	Tool for 1S163 expulsion Valve	1S164

[CONTACT US](#)
[ORDER](#)

LOCKS & VALVES

SOCKETS

Valves for Thermoformable Transfemoral Socket

	Use	Item N°.	Valve
	Valve Seat for Thermoformed Socket This valve seat will be placed on the positive model and fixed while socket thermoforming The assembly is sealed by a joint	1X21	
	Removable Valve Seat This seat valve will be fitted on thermoformable rigid plastics or soft plastics which will be laminated afterwards Use 1X4750 tool to set up	1X47	 1X16: Ø 30mm Silicone Valve
	Valve Seat and Removable Seat This seat valve will be fitted on thermoformable rigid plastics or soft plastics which will be laminated afterwards Use 1X10150 tool to set up 1X10158 wrench allows to screw the valve to the seat valve 1X10159 wrench allows to fit the seat valve on the socket	1X10140	Included

Accessories for Transfemoral Valves

	Use
	1X45 Rigid Hole Shuttle Housing To be glued on a wooden or laminated traditional socket

[CONTACT US](#)
[ORDER](#)

LOCKS & VALVES

SOCKETS

Locks for Distal Attachment

	Use	Material	Weight	Max. Patient Weight	Item N°.
	Push Lock Kit Our best selling puck-style lock with a 2.50" pin Includes (2) 2.5" lock pins and (4) M6x1.0x25 water resistant. FHCS	Plastic + Aluminum	102g	365lbs 166kg	LOC-APC-2.50
	Push Lock Kit Our best selling puck-style lock with a 1.75" pin Includes (2) 1.75" lock pins and (4) M6x1.0x25 water resistant. FHCS	Plastic + Aluminum	102g	365lbs 166kg	LOC-APC-1.75
	Push Lock Kit Our best selling puck-style lock with 1.25" pin Includes (2) 1.725" lock pins and (4) M6x1.0x25 water resistant. FHCS	Plastic + Aluminum	102g	365lbs 166kg	LOC-APC-1.25

[CONTACT US](#)
[ORDER](#)

LOCKS & VALVES

SOCKETS

Locks for Distal Attachment

	Use	Material	Weight	Tools & Accessories		Item N°.
	Locking Kit with Smooth Pin, Small Cup and Adapter for 1K90 This kit with push button and smooth pin is fitted with a device that locks the pin, allowing accurate positioning. Thanks to its small cup (Ø 57mm), it is particularly appropriate for residual limbs with low volume, such as transtibial ones. It is designed for assembly with 1K90 cup which is not included in the kit.	Plastic	100g	1S129: Tools for Lamination / Thermoforming 1S112: Centering Tool		1S137-1
	Locking Kit with Smooth Pin and Vertical Pulling This locking kit with pin is operated by vertical tab pulling. Control of unlocking can be displaced through a handle fitted on the side of the socket. 2 pins are supplied: A smooth pin that provides fine adjustment of the liner, and a smooth pin with notched end allowing the patient to know if the pin is well engaged and providing accurate adjustment. Ø 45mm cup.	Plastic + Aluminum	85g	1S141: Tools for Lamination / Thermoforming 1S129: Tool for Lamination / Thermoforming, used with 1S142 1S112: Centering Tool 1S142: Adapting Kit for assembly with 1K90 cup 1X110: Handle for displaced control of unlocking		1S140
	Locking Kit with Notched Pin and Rack Adjustment 1S145 locking kit provides very fine adjustment thanks to its rack system which avoids any slack in the bottom of the socket. It comes with a lamination prong for the upper connection and the tools needed for lamination. The lower connection is provided by screwing a threaded pyramid insert (1K182) or a threaded Receiver insert (1K183).	Aluminum + Stainless steel	125g	1K182: Threaded Pyramid insert 1K183: Threaded Receiver insert		1S145
	Smooth Pin	-	-	Pin length	30mm	1S11632
					36mm	1S11636-1
					50mm	1S11633



CONTACT US

ORDER

LOCKS & VALVES

SOCKETS

Tools and Accessories for Locks

	Use	Material	Weight	Item N°.
	European 4-Hole plate This cup allows to fit any European 4-hole linking component (4 x M6, 36mm) on a knee. Comes with its protective tools.	Aluminum	75g	1S105
	Centering Tool Fitted on the end of a positive model previously slightly leveled, the 1S112 tool allows to install a pin lock. It will be taken off (for re-use) when the positive model is broken down. Used for item No 1S103, 1S116-1, 1S117-1, 1S136-1, 1S137-1, 1S140	Plastic	-	1S112
	Forming Tool for 1S116-1 1S118-1 tool allows to keep enough space at the bottom of the socket to receive 1S116-1 locking kit with smooth pin.	Plastic	-	1S118-1
	Forming Tool for 1S117-1 1S119-1 tool allows to keep enough space at the bottom of the socket to receive 1S117-1 locking kit with smooth pin.	Plastic	-	1S119-1
	Forming Tool for 1S136-1 1S128 tool allows to keep enough space at the bottom of the socket to receive 1S136-1 locking kit with smooth pin and adapter.	Plastic	-	1S128
	Tools for 1S140 This tooling kit allows to keep enough space at the bottom of the socket to receive 1S140 locking kit. It comes with a short and a long pipe used according to the type of endoskeletal or exoskeletal socket. A drilling and screwing set is also supplied for easier assembling of the different items.	Special treated steel		1S141



CONTACT US

ORDER

LOCKS & VALVES

SOCKETS

Tools and Accessories for Locks

	Use	Material	Weight	Item N°.
	Adapting Kit for Assembly on a Cup IS142 adapting kit allows to fit IS140 locking kit with pin on 1K90 cup (not included in the kit). This kit includes a Ø 57mm cup and the screws needed to assemble the items. The use of IS129 tool is necessary.	Plastic	15g	1S142
	Socket Mounted Mechanical Knee Lock Release 1X110 handle allows to displace control of unlocking to a place that is more accessible for the patient. It is appropriate for 1M102V and 1M102V-HD locking knees.	Plastic	25g	1X110
	Threaded Pyramid Insert This pyramid is screwed to the lower end of 1S145 lock.	Plastic	65g	1K182
	Threaded Receiver Insert for 1S145 Lock This receiver is screwed to the lower end of 1S145 lock. Threaded Receiver Insert for 1S145 Lock	Plastic	75g	1X183



CONTACT US

ORDER

COVERS

SOCKETS

Transfemoral 2-part Cosmetic Cover

- Does not interfere with knee unit or its settings
- Lets clothing slip on the cosmetic cover without snagging it



TECHNICAL FEATURES

Shin Item N°	Knee Cap Item No	Knee
1G21	-	1M01 / 1M10 / 1M01-HD / 1M10-HD
1G21	1M10294	1M102 / 1M102V / 1M102-HD / 1M102V-HD

JUNCTION PLATES



Item N°	Foot
1D20	S.A.C.H.
1D26	Child



COVERS

SOCKETS

Polyethylene & Polyurethane One-part Cosmetic Covers

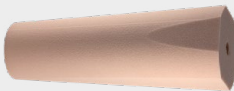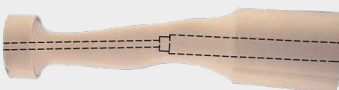	Item N°		Length	Color	Material	Diameter	Prosthesis	Packaging
	1G14		545mm	Light	Polyurethane	max 200mm	BK	x 1
	1G17		480mm		Polyethylene	Shin 147mm Ankle 100mm	BK	x 1
	1G17#12							x 12
	Left	1G13-G32	900mm		Polyurethane	320mm	1M01 1M01-HD 1M10 1M10-HD 1M102 1M102-HD 1M102V 1M102V-HD 1M03 1M05 1M112 1M113 1M115 1M149	x 1
	Right	1G13-D32						
	Left	1G13-G36				360mm		
	Right	1G13-D36						
	Left	1G13-G40				400mm		
	Right	1G13-D40						
	Left	1G13-G44				440mm		
Right	1G13-D44							
	Left	1G18-G36			Polyurethane	360mm	1P200 1P350 1P130	x 1
	Right	1G18-D36						
	Left	1G18-G40				400mm		
	Right	1G18-D40						
	Left	1G18-G44						
	Right	1G18-D44			440mm			

[CONTACT US](#)
[ORDER](#)

COVERS

SOCKETS

Plastazote One-part Cosmetic Covers

Cosmetic covers are thermoformable at 100°.

	Prosthesis	Length	Color	Material	Diameter	Item N°	Packaging
	BK	590mm	Light	Plastazote	Proximal 170mm Distal 100mm	1G1510	x 1
	AK	245mm	Medium	Polyurethane	Max 235mm	1G1525	x 1
	AK	510mm	Light	Plastazote	Proximal 190mm Distal 125mm	1G24	x 1

[CONTACT US](#)
[ORDER](#)



Materials

**Every movement depends on what supports it.
Strength where strength is needed. Flexibility
where flexibility matters. Reliability over time.**

Materials balance strength, flexibility, and durability with intention.
Weight is reduced where possible, resilience is built where needed,
and performance remains consistent over time. Every element
contributes to motion that can be trusted.

MATERIALS

TABLE OF CONTENTS



ANKLE JOINTS



STRAPS AND PADS

- Textured and colored straps and pads
- Customized pvc strap
- Fabric strap
- Velvet strap
- Micro-hook strap
- Hook and loop strap



FABRICS

- Special thermoformable fabric
- Thermoformable thick fabric
- 3D fabric
- Single layer 3d fabric waterproof
- Ripstop suede fabric
- Neoprene sheet
- Anti-slipping lining
- Waterproof & stretchable fabric
- Lycra® fabric
- Terry cloth
- Self-adhesive suede fabric
- Self-adhesive fabric
- Self-adhesive carbon imitation
- Padding felt
- Elastic velvet edge
- Thread



FASTENERS



TOOLS

[CONTACT US](#)

[ORDER](#)

ANKLE JOINTS

MATERIALS



A complete line of flexure joints, let Aflex help fabricate your all your thermoformed orthoses

- Indications: to treat pathologies of neurological, muscular or traumatic origin that affect foot levator muscles
- Easy implementation thanks to their flexibility allowing self-alignment
- Designed to last, Aflex has been tested in all conditions and guaranteed meet your needs

Flexible AFLEX

Recommended for night AFOs, their flexibility ensures comfort for even very spastic patients and possible motions of the ankle allow a better compliance of the orthosis.
(3 sizes available (pediatric, child, adult))



Description	Pediatric	Child	Adult
L x W	34 x 11,6mm	42 x 11,7mm	49 x 14mm
Thickness	5mm	6,5mm	8mm
Item N°. kit	2C131-P	2C131	2C130

Reinforced Straight AFLEX


Recommended for walking AFOs, they prevent plantarflexion of the foot and motion in rotation and inversion-eversion of the foot. Reinforced with a Kevlar fiber braid: Superior Tensile Strength (2 sizes available child, adult).



Description	Child	Adult
L x W	40 x 14mm	44 x 16mm
Thickness	9mm	
Item N°. kit	2C161	2C160

Reinforced Dorsi AFLEX

Recommended for walking AFOs with dorsiflexion assist. The natural yield of the joint provides a motion assistance in dorsiflexion during swing phase to assist toe clearance. Reinforced with a Kevlar fiber braid: Superior Tensile Strength. 3 durometers available to give a better answer for patient needs (70, 80 and 90 ShA). 2 sizes available (child, adult).

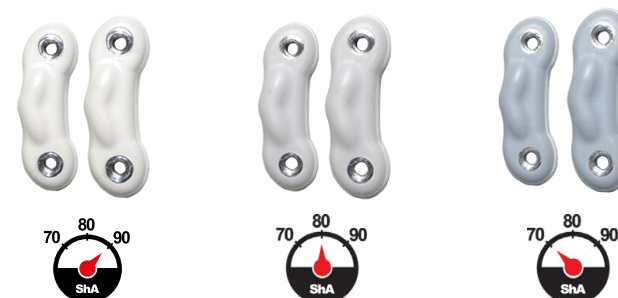

Description	White 90 ShA		Light grey 80 ShA		Dark grey 70 ShA	
	Child	Adult	Child	Adult	Child	Adult
L x W	40 x 14mm	44 x 16mm	40 x 14mm	44 x 16mm	40 x 14mm	44 x 16mm
Thickness	9mm					
Item N°. kit	2C163-90	2C162-90	2C163-80	2C162-80	2C163-70	2C162-70



CONTACT US

ORDER

STRAPS AND PADS

MATERIALS

A Full Range of Straps and Pads, Including Chafe and Loop

		SELF-GRIPPING STRAP						PVC PAD	
Pack.	Length	300mm		500mm	400mm	680mm	150mm		
		Width	16mm	25mm	40mm		25mm	40mm	
TEXTURED	x 10	Paris	2X716-CV300	2X416-CV300	2X416-CV500	2X516-CV400	2X516-CV680	2X814	2X614
	x 10	Texas	2X79-CV300	2X49-CV300	2X49-CV500	2X59-CV400	2X59-CV680	2X89	2X69
	x 100	Texas marked	2X79-CC300	2X49-CC300	2X49-CC500	2X59-CC400	2X59-CC680		
	x 10	Sevilla	2X709-CV300	2X409-CV300	2X409-CV500	2X509-CV400	2X509-CV680	2X805	2X605
	x 100	Sevilla marked	2X709-CC300	2X409-CC300	2X409-CC500	2X509-CC400	2X509-CC680		
	x 10	Sydney	2X701-CV300	2X401-CV300	2X401-CV500	2X501-CV400	2X501-CV680	2X801	2X601
	x 10	New Dehli	2X702-CV300	2X402-CV300	2X402-CV500	2X502-CV400	2X502-CV680	2X802	2X602
	x 10	Tokyo	2X70-CV300	2X40-CV300	2X40-CV500	2X50-CV400	2X50-CV680	2X85	2X65
	x 100	Tokyo marked	2X70-CC300	2X40-CC300	2X40-CC500	2X50-CC400	2X50-CC680		
	COLORED	x 10	Alaska	2X75-CV300	2X45-CV300	2X45-CV500	2X55-CV400	2X55-CV680	2X81
x 10		Moscou	2X713-CV300	2X413-CV300	2X413-CV500	2X513-CV400	2X513-CV680	2X818	2X618
x 10		Oslo	2X719-CV300	2X419-CV300	2X419-CV500	2X519-CV400	2X519-CV680	2X83	2X63
x 100		Oslo Marked	2X719-CC300	2X419-CC300	2X419-CC500	2X519-CC400	2X519-CC680		
x 10		Shanghai	2X721-CV300	2X421-CV300	2X421-CV500	2X521-CV400	2X521-CV680	2X809	2X609
x 100		Shanghai marked	2X721-CC300	2X421-CC300	2X421-CC500	2X521-CC400	2X521-CC680		
x 10		Auckland	2X714-CV300	2X414-CV300	2X414-CV500	2X514-CV400	2X514-CV680	2X819	2X619
x 100		Auckland marked	2X714-CC300	2X414-CC300	2X414-CC500	2X514-CC400	2X514-CC680		
x 10		Amsterdam	2X708-CV300	2X408-CV300	2X408-CV500	2X508-CV400	2X508-CV680	2X808	2X608
x 100		Amsterdam marked	2X708-CC300	2X408-CC300	2X408-CC500	2X508-CC400	2X508-CC680		
	x 10	Lima	2X71-CV300	2X41-CV300	2X41-CV500	2X51-CV400	2X51-CV680	150mm 25mm 150mm 40mm	
	x 100	Lima marked	2X71-CC300	2X41-CC300	2X41-CC500	2X51-CC400	2X51-CC680		
	x 10	London	2X77-CV300	2X47-CV300	2X47-CV500	2X57-CV400	2X57-CV680		

TEXTURED

COLORED

Paris

Texas

Sevilla

Sydney

New Dehli

Tokyo

Alaska

Moscou

Oslo

Shanghai

Auckland

Amsterdam

Lima


London

CUSTOMIZATION

CUSTOMIZATION

100 pieces minimum in same size and color.

Please, contact us to get an estimate, from 100 straps.

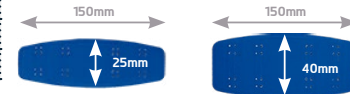

STRAPS AND PADS

MATERIALS



Customized PVC Strap

Pack.	SELF-GRIPPING STRAP			PVC PAD	
	Length	300mm	400mm	150mm	
	Width	25mm	40mm	25mm	40mm
x 10	Ladybug	2X412-CV300R	2X512-CV400R	2X812	2X612
x 10	Syxtty Flower	2X411-CV300R	2X511-CV400	2X811	2X611
x 10	Jean	2X410-CV300R	2X510-CV400R	2X810	2X610

Fabric Strap

- Flexible, comfortable and washable
- Ideal for areas sensitive to rubbing
- One size. Can be easily cut with scissors
- Delivered with roller chafe to facilitate the sliding of the strap

Pack.	Length	290 + 125mm	290 + 125mm
	Width	25mm	40mm
x 10	Item N°.	2X100-2525	2X100-4025

Velvet Strap

- Flexible, comfortable and washable strap
- It is ideal for friction sensitive areas
- Better comfort thanks to the built-in flap
- Impressive gripping strength thanks to the micro hook
- Delivered with roller chafe to facilitate the sliding of the strap

Pack.	Width	25mm	40mm
	Item N°.	2X098-2525	2X098-4025
x 10			

Micro-Hook Strap

- Velvet-hemmed loop band, associated with our micro hook pliers
- Optimal implementation: can be cut to the desired length
- Delivered with roller chafe to facilitate the sliding of the strap

Pack.	Width	25mm	40mm
	Item N°.	2X099-25	2X099-40
x 10			




CONTACT US

ORDER

FABRICS

MATERIALS

Hook and Loop Strap

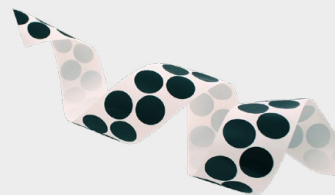							
Comfort Pad		Elastic Loop Strap		Micro-hook		Plastic Disc with Hooks	
- Provides comfort in areas sensitive to rubbing. - The hooks on one side make it easy to attach to loop. - 2X100-NS with silicone to prevent the orthosis from slipping. - Easy cutting. - Particularly adapted to 2x100-4025		- Adjustable and variable compression. - Perfectly matches the skin relief. - Easy cutting.		- Very flat, soft to the touch. - Exceptional gripping capacity. - Bonds easily to many supports like polyethylene, polypropylene, PVC, polyamide, wood, glass, cardboard, metal. - Precut round pieces. - Adhesive side for immediate bond, optimal after 24h.		- Very strong grip. - Can be screwed in the plastic device (screws included).	
Width	50mm		25mm	50mm	Ø 15mm	Ø 25mm	
Length	1m		25m		-		-
Color	Black	Black with silicone	White		Black	Black	Transparent
Item N°.	2X100-N	2X100-NS	TT190-25	TT190-50	TT216-15	TT216-25	TT215-25
Packaging	by meter	10m	x 1		x 800	x 400	25m
							10 discs + 10 screws

[CONTACT US](#)
[ORDER](#)

FABRICS

MATERIALS



Special Thermoformable Fabric

- Breathable, textured, moisture wicking fabric
- Very comfortable and soft to the touch
- Easy to clean, washable at 40°C
- Good resistance to aging
- Glue-less, this fabric limits risks of allergy
- Thermoformable at 170°C

Thickness	6mm	Width	1.45m	Length	1.00m
-----------	-----	-------	-------	--------	-------

							
Black	White	Blue	Raspberry	Lead gray	Brown	Green	Red
TT501	TT502	TT503	TT504	TT505	TT506	TT507	TT508

Thermoformable Thick Fabric

- Perfectly matches the most complex mold shapes
- Very comfortable and soft to the touch
- Easy to clean, washable at 40°C
- Good resistance to aging
- Glue-less, this fabric limits risks of allergy
- Thermoformable at 170°C

Thickness	6mm	Width	1.45m	Length	1.00 m
-----------	-----	-------	-------	--------	--------

		
Black	Blue	Fuchsia
PL190-06	PL191-06	PL192-06

CONTACT US

ORDER

A stack of various colored, perforated mats. The colors from top to bottom are black, white, blue, pink, grey, brown, green, and red. A circular inset shows a cross-section of the material, revealing a porous, fibrous structure with arrows indicating air flow or moisture wicking.

- Breathable, textured, moisture wicking fabric
- Very comfortable and soft to the touch
- Easy to clean, washable at 40°C
- Good resistance to aging
- Glue-less, this fabric limits risks of allergy

Thickness:	6mm	Width	1.45m	Length	by meter		
							
Black	White	Blue	Raspberry	Lead gray	Brown	Green	Red
TT381	TT382	TT383	TT384	TT385	TT386	TT387	TT388

- Stretchable, breathable, comfortable 3D fabric
- Waterproof thanks to the associated flexible film
- Very comfortable and soft to the touch
- Ideal for making covers with waterproof areas



Thickness	3mm	Width	1.45m	Length	by meter
					
Black	White	Blue	Raspberry	Lead gray	Red
TT391	TT392	TT393	TT394	TT395	TT396



ORDER

FABRICS

MATERIALS

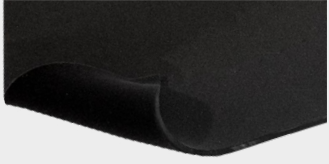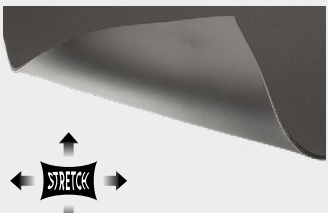									
Ripstop Suede Fabric		Neoprene Sheet		Anti-slipping Lining		Waterproof & Stretchable Fabric		Lycra® Fabric	
- Its texture allows it to have strong resistance and velvet soft touch. - Very comfortable contact, soft and flexible to the touch. - Easy to clean, washable up to 40°C. - Good resistance to aging.		- Particularly suitable for lining orthosis and the confection of plastron. - Neoprene sheet with one side in polyester fabric, the back side can be fixed with hook. - Easy to clean, washable up to 40°C.		- Particularly suitable for lining of orthosis closure bands. - This copolymer provides optimal adhesion to the patient skin, the back side can be fixed with hook. - Relieves friction spots. - Easy to clean with water and soap.		- Particularly suitable for making waterproof covers. - Extensibility (bi stretch). - Excellent mechanical strength. - Easy to clean, washable at 40°C for TT372 and up to 90°C for TT373#10.		- Particularly suitable for making covers. - Extensibility (bi stretch). - Excellent mechanical strength. - Easy to clean, washable at 40°C. - Dries faster than an ordinary fabric.	
Thickness	4,5mm		5mm		5mm		2,5mm	1mm	-
Width	1,45m		1,45m		0,30m		1,45mm	1,45mm	1,5m
Length	by meter	10m	0,99m		0,30m		by meter	10m	10m
Item N°.	TT600	TT600#10	TT460		TT450		TT372	TT373#10	TT370#10
Packaging	by meter	10m	x 1		x 1		by meter	10m	10m

[CONTACT US](#)
[ORDER](#)


FABRICS

MATERIALS




Terry cloth

- Particularly suitable for making covers.
- Ideal for lining.
- Stretchable in both directions.
- Easy to clean, washable up to 40°.

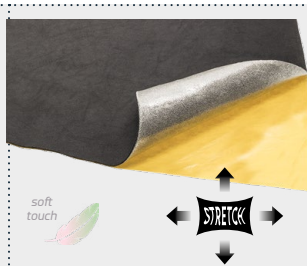

Self-adhesive Suede Fabric

- It is recommended for inside or outside lining of orthosis.
- Very comfortable contact, soft and flexible to the touch.
- Fabric with a powerful adhesive.
- Good resistance to aging.

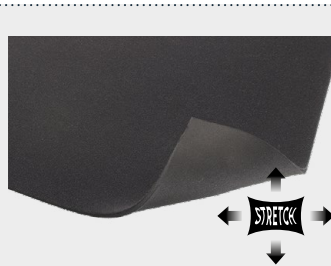

Self-adhesive Fabric

- Ideal for inside or outside lining of orthosis.
- Very comfortable and soft to the touch.
- Fabric with a powerful adhesive.



Self-adhesive Carbon Imitation

- This vinyl is particularly recommended for the lining of orthosis.
- It allows the covering of any material in order to get a «carbon-fiber» appearance.
- Resistant and easy to apply (with a heat gun).
- Adhesive side for immediate bond, optimal after 24h.

Padding Felt

- Particularly suitable for making covers.
- Extensibility (bi stretch).
- Excellent mechanical strength.
- Easy to clean, washable at 40°C.
- Dries faster than an ordinary fabric.

Thickness	2mm			1,5mm	1,5mm	-	2mm
Width	1,5m			1,5m	1,5m	1,52m	0,9m
Length	5m			by meter	by meter	1m	1m
Color	Blue ●	Pink ●	Black ●	Black ●	Black ●	Carbon imitation	White ○
Item N°.	TT550-5	TT551-5	TT552-5	TT410	TT400	TH001	TF001
Packaging	5 m			x 1	by meter	x 1	x 1

CONTACT US

ORDER

FABRICS

MATERIALS



Elastic Velvet Edge

- Its elasticity makes it suitable for finishing of edges with complex shapes
- Very comfortable and soft to the touch
- Reduces risks of irritation caused by rough edges
- Can be stitched or glued on plastics

										
Color	Black	White	Blue	Fuchsia	Grey	Brown	Red	Green	Beige	
Width	25mm								22mm	
Packaging	25m									
Item N°.	TT020-N	TT020-BLANC	TT020-BLEU	TT020-F	TT020-G	TT020-M	TT020-R	TT020-V	TT020	TT025

Thread



											
Color	White			Raw	Black	Blue	Fuchsia	Grey	Brown	Green	Red
Material	Linen	Terylene	Polyester	Polyester	Polyester						
Packaging	50g	270g	235g	250g	235g						
Item N°.	TT050	TT055	TT056	TT051-1	TT056-N	TT056-B	TT056-F	TT056-G	TT056-M	TT056-V	TT056-R

[CONTACT US](#)
[ORDER](#)

FASTENERS

MATERIALS

Colored Single Cap Speedy Rivet

Head	Ø 13mm	Tube	Ø 4mm	Pack.	x 100
------	--------	------	-------	-------	-------

																
Color	Black		White		Blue		Fuchsia		Grey		Brown		Green		Red	
Height	15,6mm	12,4mm	15,6mm	12,4mm	15,6mm	12,4mm	15,6mm	12,4mm	15,6mm	12,4mm	15,6mm	12,4mm	15,6mm	12,4mm	15,6mm	12,4mm
Item N°.	KE100	KE200	KE101	KE201	KE102	KE202	KE103	KE203	KE104	KE204	KE105	KE205	KE106	KE206	KE107	KE207

Single Cap Speedy Rivet

Packaging	Head Ø	6,6mm		8,5mm			9,5mm	10,5mm	12mm
x 100	Tube Ø	2,5mm		3mm			3mm	3mm	3mm
	Height	6mm	7mm	8mm	11mm	15,6mm	9,1mm	11mm	10mm
	Item N°.	KE047	KE046	KE050	KE051	KE057	KE056	KE058	KE053



Double Cap Speedy Rivet

Packaging	Head Ø	10,7mm
	Tube Ø	3mm
x 100	Height	9,1mm
	Item N°.	KE052



Press Studs (snap)

Packaging						
	Color	Steel	Black	White	Blue	Fuchsia
x 100	Item N°.	KE030	KE030-01	KE030-02	KE030-03	KE030-04


[CONTACT US](#)
[ORDER](#)

TOOLS

MATERIALS



Scissors

Description	Length	Item N°.
Plaster cast scissors	240mm	5C100
Scissors	700mm	5C101
Dressmaker's scissors	170mm	5C110
Dressmaker's scissors	200mm	5C111

Deburring Tool

Ideal for trimming and leveling edges and surfaces

Description	Item N°.
Trimming with 2 blades	5C425
10 spare blades	5C427

Deburring Tool for Pylons

One side to trim inside and the other to trim outside

Description	Item N°.
Inside and outside trimming for 22 to 34mm Ø tubes	5C424

Metal Forming Form

Jaw opening	Length	Item N°.
4 & 6mm	245mm	5D210
	310mm	5D220
6 & 8mm	245mm	5D211
	310mm	5D221
7 & 9mm	310mm	5D222

[CONTACT US](#)
[ORDER](#)

PROTEOR TERMS OF SALE, RETURNS, & WARRANTIES

TABLE OF CONTENTS

1.0 NEW ACCOUNTS

1.1 How to Open an Account:

We establish an open account when you place your first order with PROTEOR. We require a completed credit application prior to processing a second order. The review and approval process will normally take no longer than 24 hours from our receipt of your application. Orders must be pre-paid or shipped C.O.D. until credit has been established. Obtain a credit application by telephone at: 855.450.7300.

1.2 How to Order:

You may place orders with PROTEOR in the USA, APAC, and Latin America by: Telephone: 855.450.7300 (Toll free, USA) Email: orders@proteor.com

1.3 Order Information

When transmitting your order, please ensure that the order includes the following information.

- Part # or Item Description
- Quantity desired
- Order need by date
- "Ship to" address and phone number
- "Bill to" address (if it differs from "Ship to" address)
- Transportation mode desired (UPS Ground, 2-day, etc.)

Upon our receipt, thereof, our Customer Service Department will confirm the order and schedule it for shipment.

2.0 RETURN POLICY & PROCESS

2.1 Returns, Direct Customers

2.1.2 Return Authorization Number (RA#) Required

A. All returns require a Return Authorization Number (RA#). Call PROTEOR's USA Customer Service Department at 855.450.7300 for authorization. Be prepared to provide the Serial Number of product(s) you seek to return.

B. Place the provided Return Authorization Number inside and on the shipping container, on or near the shipping label. We will not accept returns without a RA#.

2.1.3 Return Destination: Authorized returns are to be shipped to the location address listed on the Return Authorization for the specific product.

2.1.4 Return Opportunity: We will not accept returns after 60 days from the invoice date. We do not accept returns for Knees or Microprocessor products.

2.1.5 Return Documentation: Send a copy of the original invoice or packing list and the customer return card with the returned product to ensure that you receive the replacement or credit for which you qualify. Failure to provide the required documentation may delay or prevent return processing.

2.1.6 Responsibility: Until PROTEOR receives the returned product, the Prosthetist remains responsible for the product and all related billing.

2.1.7 Custom Orders: Prosthetist is responsible for all shipping costs. We do not accept returns for orders (custom products). Custom orders are products of a Line, Style, Category, or Size not ordinarily stocked in our distribution facility.

2.1.8 Products That Qualify for Return between 0 - 60 Days of Invoice Date:

- | | |
|-------------------------|------------------------------|
| • Kinterra | • EasyRun |
| • ShockWave | • Gery |
| • RUSH ROGUE 2 | • DynaStep |
| • RUSH HiPro | • Foot Shells |
| • RUSH Rover | • Components |
| • Highlander | • Liners and Sealing Sleeves |
| • RUSH Chopart | • Sockets |
| • Agilix | • Locks and Valves |
| • ETHOS LP | • Cosmetics |
| • Pacifica LP | • Ankle Joints |
| • DynAdapt | • Straps and Pads |
| • RUSH Kid | • Fabrics |
| • RUSH H2O Collection | • Fasteners |
| • RUSH EVAQ8 Collection | • Tools |

2.1.9 Resaleable Condition: All products returned within 60 days of the invoice date will be evaluated to determine that they are in resaleable condition. Such products found to be damaged or not in resaleable condition shall be returned to the Customer, with no credit issued.

PROTEOR TERMS OF SALE, RETURNS, & WARRANTIES

TABLE OF CONTENTS

2.2 Types of Returns

2.2.1 60-day Returns: If a Customer is dissatisfied with standard stock products that are eligible for a Trial Period, they may return the product for size alteration, stiffness modification, or account credit if product is in resalable condition.

A. Exceptions:

- Special order sizes, category variations, and custom foot products are not eligible for returns.
- Microprocessor products are non-returnable: QUATTRO, Plie, SYNSYS, and Kinnex
- Knees are non-returnable.

2.2.2 Miscellaneous “No Replacement” Returns: When a Customer does not place a reorder at the time the RA# is requested.

- We will accept all returns with RA#s issued between 0 and 60 days of the invoice date and will not charge a restocking fee if the product is in resaleable condition.

2.3 Account Credit: We grant an account credit and apply it to the Customer’s account, pursuant to existing procedure or a negotiated agreement with the Customer (e.g., credit for returned goods, overpayments, duplicate payments, etc.).

3.0 WARRANTY POLICY & PROCEDURES

3.1 Warranty Policy: We guarantee all our products to be free from defects in material and workmanship. Any product found to be defective, when used in a customary and proper manner according to the product’s published Instructions for Use, shall at PROTEOR’s discretion be evaluated and replaced under warranty. All warranty replacements will be shipped in a timely manner accompanied by an invoice for the original purchase price. We will issue credit upon receipt of the product’s return to PROTEOR.

3.2 Warranty Terms: We calculate Warranty duration from the original invoice date.

3.2.1 Ankle / Foot Modules: 36-months

A. Exceptions:

- | | |
|------------------------|--------------------------------------|
| • Gery Foot: 12-months | • ETHOS LP for Pediatrics: 18-months |
| • DynaStep: 12-months | • RUSH Chopart: 24-months |
| • EasyRun: 12-months | • Hopper: 24-months |

3.2.2 Foot Shells: 6-months

3.2.3 Microprocessor Ankles

- **PROTEOR Kinnex:** No required service checks. Extended 2- year warranty option available when purchased within 12 months of the original invoice date.

3.2.4 Microprocessor Knees: 36-months if mandatory service checks are met.

- **PROTEOR QUATTRO:** No required service checks. Extended 2- and 3-year warranty options available when purchased within 12 months of the original invoice date.
- **Plie:** Mandatory service checks at 12 and 24 months. Extended 2-year warranty is available when purchased within 12 months of the original invoice date.
- **SYNSYS:** 36 months warranty with mandatory service at 24 months. Extended 3 year warranty when purchased within 12 months of the original invoice date and mandatory servicing at 48 months.

3.2.5 Mechanical Knees: 36-months

A. Exceptions:

- Symphony: 24-months
- EasyRide: 12-months

3.2.6 Structural Components: 36-months

3.2.7 Socket Technology: No Warranty

A. Exception:

- Keasy Cones: 24-months
- SiLiTouch Liners: 6-months

3.2.8 Orthotics:

- Straps: 6-months
- Joints: 12-months

3.2.9 Fabrication Supplies: 3-months

3.3 Exclusions: The warranties set forth herein do not cover the battery, battery charger, adapters and optional accessories (each of which are covered under separate warranty), or test sockets, nor do they cover any alteration, repair or substitution of other parts at PROTEOR judgment, materially or adversely affects the product or part, damage due to accidents, cosmetic damage, water damage deemed as abuse, neglect, misuse or operation beyond capacity, or parts damaged by improper installation.

DISCLOSURE: THE FOREGOING WARRANTIES ARE IN LIEU OF ALL WARRANTIES, EITHER EXPRESS OR IMPLIED, WRITTEN OR ORAL, INCLUDING WITHOUT LIMITATION, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, OR OTHERWISE, AND OF ANY OTHER OBLIGATION ON THE PART OF SUPPLIER. THE SOLE AND EXCLUSIVE LIABILITY OF SUPPLIER HEREUNDER SHALL BE TO REPLACE, REPAIR OR REFUND THE PURCHASE PRICE OF A DEFECTIVE PRODUCT, AS DESCRIBED ABOVE. NO WAIVER, ALTERATION, OR MODIFICATION OF THE FOREGOING WARRANTIES SHALL BE BINDING AGAINST SUPPLIER, UNLESS IT IS IN WRITING AND SIGNED BY AN AUTHORIZED EXECUTIVE OFFICER OF SUPPLIER.

Terms of Sales, Returns, and Warranties – February 2025

TERMS OF USE

TABLE OF CONTENTS

EFFECTIVE DATE: JANUARY 1, 2021

ACCEPTANCE OF TERMS OF USE

These terms of use are entered into by and between you and Proteor USA, LLC ("Company", "we", "our" or "us"). The following terms and conditions, together with any documents they expressly incorporate by reference (collectively, these "Terms" or "Agreement"), govern your access to and use of the Company's website, including any content, functionality, and services offered on or through us.proteor.com (the "Site").

Please read these Terms carefully before you start to use the Site. By (i) accessing, browsing, subscribing to, and using the Site; and (ii) by using our associated content and any customizations made to services provided or offered by us via the Site (collectively, the "Services"), you acknowledge that you have read, have understood, and you accept and agree to be bound by this Agreement, to abide by these Terms, and to comply with all applicable laws and regulations. If you do not agree to these terms, you must not access or use the Site or any of the Services.

These Terms apply to you and your authorized users. If you accept these Terms on behalf of an entity, you hereby represent and warrant that you are authorized to do so. You may permit your employees, agents, or contractors ("Authorized Users") to use the Services, provided that your authorized users comply with these Terms. You are responsible for the use and actions of your authorized users.

This Site is offered and available to users who are 18 years of age or older. By using this Site, you represent and warrant that you are of legal age to form a binding contract with us and meet all of the foregoing eligibility requirements. If you do not meet all of these requirements, you must not access or use the Site.

The material provided through our Site and the Services is protected by law, including, but not limited to, United States copyright law and international treaties. We operate the Site and offer our Services from our offices within the United States. We make no representation that materials on our Site and in our Services are appropriate or available for use in other locations, and access to them from territories where their contents are illegal is prohibited. Those who choose to access the Site and the Services from other locations do so on their own initiative and are responsible for compliance with all applicable local laws.

CHANGES TO THE TERMS

Information provided on the Site and the Services may be changed or updated without notice. Company may also make improvements and/or changes in the Services described in this information at any time without notice. Your continued use of the Site following the posting of revised Terms means that you accept and agree to the changes. You are expected to check this page from time to time so you are aware of any changes, as they are binding on you.

ACCESSING THE SITE, REGISTRATION, ACCOUNTS, AND PASSWORDS

We reserve the right to withdraw or amend the Site, and any service or material we provide on the Site, in our sole discretion without notice. We will not be liable if for any reason all or any part of the Site is unavailable at any time or for any period. From time to time, we may restrict access to some parts

of the Site, or the entire Site, to users, including registered users. You are responsible for making all arrangements necessary for you to have access to the Site and ensuring that all persons who access the Site through your internet connection are aware of these Terms and comply with them.

USE RESTRICTIONS AND INTELLECTUAL PROPERTY

The Site and its entire contents, features, and functionality (including but not limited to all information, software, text, displays, images, video, and audio, and the design, selection, and arrangement thereof ("Company Content") is owned by the Company, its licensors, or other providers of such material (collectively, and individually, "Company Content Providers") and are protected by United States and international copyright, trademark, patent, trade secret, and other intellectual property or proprietary rights laws. You agree not to modify, copy, distribute, transmit, share, display, perform, reproduce, publish, license, create derivative works from, transfer, or sell or re-sell any Company Content obtained from or through Company. The Site and Services contain proprietary and copyright-protected information. Any unauthorized use of any Company Content contained on the Site or the Services may violate copyright laws, trademark laws, the laws of privacy and publicity, and communications regulations and statutes. You may not modify, publish, transmit, participate in the transfer or sale, create derivative works, or in any way exploit any of the content in whole or in part.

You will not upload, post or otherwise make available on the Site or via the Services, any material protected by copyright, trademark or other proprietary right without the express permission of the owner of the copyright, trademark or other proprietary right. You will indemnify Company and its officers, directors, employees, and agents for any claims by third parties of infringement of copyrights, trademarks, proprietary rights or any other harm resulting from such a submission of material by you.

The Company Content and Site may not be reproduced, duplicated, copied, sold, resold, or otherwise exploited for any commercial purpose without express written consent of Company. You may not frame or utilize framing techniques to enclose any trademark, logo, or other proprietary information (including images, text, page layout, or form) of Company without express written consent. You must not modify copies of any materials from this site or delete or alter any copyright, trademark, or other proprietary rights notices from copies of materials from this site. You may not use any meta tags or any other "hidden text" utilizing Company's name or trademarks without the express written consent of Company. You may not misuse the Site. You may use the Site only as permitted by law. You must not access or use for any commercial purposes any part of the Site or any services or materials available through the Site.

The Company name, the term PROTEOR, PROTEOR USA, the Company logo, and all related names, logos, product and service names, designs, and slogans are trademarks of the Company or its affiliates or licensors. You must not use such marks without the prior written permission of the Company. All other names, logos, product and service names, designs, and slogans on this Site are the trademarks of their respective owners.

No right, title, or interest in or to the Site or any content on the Site is transferred to you, and all rights not expressly granted are reserved by the Company. Any use of the Site not expressly permitted by these Terms is a breach of these Terms and may violate copyright, trademark, and other laws.

TERMS OF USE

TABLE OF CONTENTS

DISCLAIMER OF WARRANTIES

THE INFORMATION PRESENTED ON OR THROUGH THE SITE IS MADE AVAILABLE SOLELY FOR GENERAL INFORMATION PURPOSES. WE DO NOT WARRANT THE ACCURACY, COMPLETENESS, OR USEFULNESS OF THIS INFORMATION. ANY RELIANCE YOU PLACE ON SUCH INFORMATION IS STRICTLY AT YOUR OWN RISK. WE DISCLAIM ALL LIABILITY AND RESPONSIBILITY ARISING FROM ANY RELIANCE PLACED ON SUCH MATERIALS BY YOU OR ANY OTHER VISITOR TO THE SITE, OR BY ANYONE WHO MAY BE INFORMED OF ANY OF ITS CONTENTS. INFORMATION ON THE SITE MAY CONTAIN TECHNICAL INACCURACIES OR TYPOGRAPHICAL ERRORS. ALL INFORMATION IS PROVIDED BY COMPANY ON AN "AS IS" BASIS ONLY. COMPANY PROVIDES NO REPRESENTATIONS AND WARRANTIES, EXPRESS OR IMPLIED, INCLUDING THE IMPLIED WARRANTIES OF FITNESS FOR A PARTICULAR PURPOSE, MERCHANTABILITY AND NON-INFRINGEMENT, UNLESS SUCH REPRESENTATIONS AND WARRANTIES ARE NOT LEGALLY EXCLUDABLE. USE OF THE SERVICES IS SOLELY AT YOUR RISK. COMPANY MAKES NO WARRANTY THAT (i) THE SITE WILL MEET YOUR REQUIREMENTS; (ii) THE SITE WILL BE UNINTERRUPTED, TIMELY, SECURE, OR ERROR FREE; (iii) THE RESULTS THAT MAY BE OBTAINED FROM THE USE OF THE SITE WILL BE ACCURATE OR RELIABLE; AND (iv) THE QUALITY OF ANY PRODUCTS, INFORMATION, OR OTHER SERVICES PURCHASED OR OBTAINED BY YOU THROUGH THE SITE WILL MEET YOUR EXPECTATIONS.

YOU ARE SOLELY RESPONSIBLE FOR ENSURING THAT YOUR USE OF THE SITE, RELATED SERVICES, OR CONTENT IS IN ACCORDANCE WITH APPLICABLE LAW. COMPANY IS NOT RESPONSIBLE FOR ANY INFORMATION NOT PROVIDED ON THE SITE. COMPANY DOES NOT GUARANTEE, IMPLICITLY OR EXPLICITLY, ANY RESULTS. ANY MATERIAL DOWNLOADED OR OTHERWISE OBTAINED THROUGH THE USE OF THE SITE IS DONE AT YOUR OWN DISCRETION AND RISK AND YOU WILL BE SOLELY RESPONSIBLE FOR ANY DAMAGE TO YOUR COMPUTER SYSTEM OR LOSS OF DATA THAT RESULTS FROM THE DOWNLOAD OF ANY SUCH MATERIAL. NO ADVICE OR INFORMATION, WHETHER ORAL OR WRITTEN, OBTAINED BY YOU FROM COMPANY OR THROUGH OR FROM THE SITE SHALL CREATE ANY WARRANTY NOT EXPRESSLY STATED IN THE TERMS OF SERVICE AND USE.

LIMITATION OF LIABILITY

YOU EXPRESSLY UNDERSTAND AND AGREE THAT COMPANY SHALL NOT BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, CONSEQUENTIAL OR EXEMPLARY DAMAGES, INCLUDING BUT NOT LIMITED TO, DAMAGES FOR LOSS OF PROFITS, REVENUES, GOODWILL, USE, DATA, OR OTHER INTANGIBLE LOSSES (EVEN IF COMPANY HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES), RESULTING FROM: (i) THE USE OR THE INABILITY TO USE THE SITE; (ii) THE COST OF PROCUREMENT OF SUBSTITUTE GOODS AND SERVICE RESULTING FROM ANY GOODS, DATA, INFORMATION, OR SERVICE PURCHASED OR OBTAINED OR TRANSACTIONS ENTERED INTO THROUGH OR FROM THE SITE; (iii) UNAUTHORIZED ACCESS TO OR ALTERATION OF YOUR TRANSMISSIONS OR DATA; (iv) STATEMENTS OR CONDUCT OF ANY THIRD PARTY ON THE SITE; OR (v) ANY OTHER MATTER RELATING TO THE SERVICES. YOU ALSO AGREE THAT COMPANY WILL NOT BE LIABLE FOR ANY INTERRUPTION OF BUSINESS, ACCESS DELAYS, OR ACCESS INTERRUPTIONS TO THE SITE OR SERVICES, DATA NON-DELIVERY, MIS-DELIVERY, CORRUPTION, DESTRUCTION OR OTHER MODIFICATION, AND EVENTS BEYOND COMPANY'S REASONABLE CONTROL.

IN NO EVENT SHALL COMPANY BE FINANCIALLY LIABLE FOR ANY CONSEQUENTIAL OR INCIDENTAL DAMAGES. BECAUSE SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF LIABILITY FOR CONSEQUENTIAL OR INCIDENTAL DAMAGES, IN SUCH STATES, OUR LIABILITY IS LIMITED TO THE EXTENT PERMITTED BY LAW.

TRANSMISSIONS AND CUSTOMER DATA

"Customer Data" means all data, information, or content provided, uploaded, recorded, or posted by any Authorized User during and in connection with its authorized use of the Services published or displayed on or through the Services or transmitted to or shared with other users. Customer Data includes protected health information (PHI) but does not include Aggregated Data. Any Customer Data, material, information, or idea you transmit to the Site or through the Services by any means will be treated in accordance with Company's Privacy Policy.

You grant us the limited license to process the Customer Data you transmit on our Site or through our Services. You represent and warrant that you own or control all rights in, to, and under the Customer Data you transmit on our Site or through our Services, and have received the applicable consent to collect, use, and process such Customer Data, and have the right to grant the license granted above to us and our affiliates and service providers, and each of their and our licensees, successors, and assigns.

You understand and acknowledge that you are responsible for any Customer Data you submit or transmit, and you, not the Company, have full responsibility for such content, including its legality, reliability, accuracy, and appropriateness.

AGGREGATED DATA

You authorize Company to compile Aggregated Data and such Aggregated Data shall be proprietary to Company. To the extent necessary, you hereby assign to Company all right, title, and interest in, to, and under the Aggregated Data.

"Aggregated Data" means any non-personally identifiable, technical, statistical or analytical data, gathered or generated directly by the Services or Deliverables or by use of the Services or Deliverables, that Company collects, gathers, and aggregates periodically as part of its ordinary business operations. Aggregated Data is de-identified data that Company (and its affiliates, licensors, and agents) may use on a non-attributed basis to monitor and improve its products and services, for benchmarking purposes, or to provide customized services or technologies to its customers

PROHIBITED USES

You may use the Site only for lawful purposes and in accordance with these Terms. You agree not to use the Site:

- In any way that violates any applicable federal, state, local, or international law or regulation (including, without limitation, any laws regarding the export of data or software to and from the US or other countries).
- For the purpose of exploiting, harming, or attempting to exploit or harm minors in any way by

TERMS OF USE

TABLE OF CONTENTS

exposing them to inappropriate content, asking for personally identifiable information, or otherwise.

- To transmit, or procure the sending of, any advertising or promotional material, including any “junk mail”, “chain letter”, “spam”, or any other similar solicitation.
- To post or transmit to or from the Site any unlawful, threatening, libelous, defamatory, obscene, scandalous, inflammatory, pornographic, or profane material, or any other material that could give rise to any civil or criminal liability.
- To impersonate or attempt to impersonate the Company, a Company employee, another user, or any other person or entity (including, without limitation, by using email addresses associated with any of the foregoing).
- To engage in any other conduct that restricts or inhibits anyone’s use or enjoyment of the Site, or which, as determined by us, may harm the Company or users of the Site or expose them to liability.
- Use the Site in any manner that could disable, overburden, damage, or impair the site or interfere with any other party’s use of the Site, including their ability to engage in real time activities through the Site.
- Use any robot, spider, or other automatic device, process, or means to access the Site for any purpose, including monitoring or copying any of the material on the Site or reverse engineer the Site or the Services.
- Use any manual process to monitor or copy any of the material on the Site or for any other unauthorized purpose without our prior written consent.
- Use any device, software, or routine that interferes with the proper working of the Site.
- Introduce any viruses, trojan horses, worms, logic bombs, or other material that is malicious or technologically harmful.
- Attempt to gain unauthorized access to, interfere with, damage, or disrupt any parts of the Site, the server on which the Site is stored, or any server, computer, or database connected to the Site.
- Attack the Site via a denial-of-service attack or a distributed denial-of-service attack.
- Otherwise attempt to interfere with the proper working of the Site.

TERMINATION

In its sole and absolute discretion, with or without notice to you, Company may suspend or terminate your use of and access to the Services, terminate your account and/or remove and discard anything transmitted by you, or information stored, sent, or received via the Services without prior notice and

for any reason, including, but not limited to: (i) any unauthorized access or use of the Services, (ii) any violation of this Agreement, or (iii) tampering with or alteration of any of the software and/or data files contained in, or accessed through, the Services. Company shall not be liable to you or any third party for any claims or damages arising out of any termination or suspension of the Site. Termination, suspension, or cancellation of the Services or your access rights to the Site shall not affect any right or relief to which Company may be entitled, at law or in equity, and all rights granted to you will automatically terminate and immediately revert to Company.

INDEMNIFICATION

You agree to defend, indemnify, and hold harmless the Company, its affiliates, licensors, and service providers, and its and their respective officers, directors, employees, contractors, agents, licensors, suppliers, successors, and assigns from and against any claims, liabilities, damages, judgments, awards, losses, costs, expenses, or fees (including reasonable attorneys’ fees) arising out of or relating to your violation of these Terms or your use of the Site, including, but not limited to, your use of the Site’s content, services, and products (other than as expressly authorized in these Terms) or your use of any information obtained from the Site.

GOVERNING LAW AND DISPUTE RESOLUTION

This Agreement and performance hereunder shall be exclusively governed by, and construed in accordance with, the laws of the state of Delaware and in accordance with the provision defined in the applicable Masters Service Agreement

ANY CAUSE OF ACTION OR CLAIM YOU MAY HAVE ARISING OUT OF OR RELATING TO THESE TERMS OF USE OR THE SITE MUST BE COMMENCED WITHIN ONE (1) YEAR AFTER THE CAUSE OF ACTION ACCRUES, OTHERWISE, SUCH CAUSE OF ACTION OR CLAIM IS PERMANENTLY BARRED.

WAIVER AND SEVERABILITY

No waiver of by the Company of any term or condition set out in these Terms shall be deemed a further or continuing waiver of such term or condition or a waiver of any other term or condition, and any failure of the Company to assert a right or provision under these Terms shall not constitute a waiver of such right or provision.

If any provision of these Terms is held by a court or other tribunal of competent jurisdiction to be invalid, illegal or unenforceable for any reason, such provision shall be eliminated or limited to the minimum extent such that the remaining provisions of the Terms will continue in full force and effect.



CONTACT INFORMATION

Customer Service

PHONE: 855.450.7300

EMAIL: info@proteor.com

HOURS: Monday-Friday, 6AM-5PM PST

Technical and Clinical Support

PHONE: 855.450.7300

EMAIL: support@proteor.com

Placing Orders

CALL: 855.450.7300

EMAIL: orders@proteor.com

ONLINE: shop.proteorusa.com

PROTEOR

Arizona - US Headquarters

1236 West Southern Avenue
Suite 101
Tempe, AZ 85282

California

46 Discovery
Suite 101
Irvine, CA 92618

Utah

491 East 350 North
Gunnison, UT 84634

Pennsylvania

5988 Coplay Road
Whitehall, PA 18052